



867 pH Module mit Touch Control

2.867.0110

High-end pH/Ionenmeter basierend auf dem 867 pH Module, inklusive 900 Touch Control, 854 iConnect und intelligenter pH-Glaselektrode «iUnitrode».

Neben Messung von pH, Temperatur, mV, Ipol, Upol und Konzentration kann das pH Module Standardadditionen (manual, dos, autos) und Liquid Handling (add, prep, empty) durchführen. Es ermöglicht die Verwendung von sowohl herkömmlichen als auch intelligenten Sensoren zum Messen. In die Software integriert ist auch ein automatischer GLP-konformer pH-Elektrodentest.

Das pH Module besitzt 2 USB-Schnittstellen zum Anschluss von Druckern, Barcodelesern oder Probenwechslern und 4 MSB-Schnittstellen für Rührer oder Dosinos (für die Zugabe von Hilfslösungen oder für die Standardaddition).

In Kombination mit dem Touch Control erfüllt es die Anforderungen von GLP und FDA 21 CFR part 11.

Lieferumfang 2.867.0110

Qt.	Order no.	Beschreibung
-----	-----------	--------------

1 PCS

1.854.0010

854 iConnect

854 iConnect - electrode cable and measuring amplifier for intelligent electrodes «iTrodes».



1 PCS

1.867.0010

867 pH Module

Module for pH/ion measurement as supplement to a Titrando or "stand-alone" in combination with a 900 Touch Control.

In addition to measurements of pH, temperature, mV, Ipol, Upol and concentration, the pH Module can perform standard additions (manual, dos, autos) and Liquid Handling (add, prep, empty). It enables the use of both conventional and intelligent sensors for measuring. Also integrated in the software is an automatic GLP-compliant pH electrode test.

The pH Module has two USB interfaces for connecting printers, barcode readers or sample changers and four MSB interfaces for stirrers or Dosinos (for the addition of auxiliary solutions or for standard addition). Both in conjunction with the Touch Control and also integrated in *tiamo*™ full (from 2.0), it is in compliance with GLP and FDA 21 CFR part 11 requirements.




1 PCS

1.900.0010

900 Touch Control

Operating unit for the Titrandos, USB Sample Processors, 856 Conductivity Module, 867 pH Module and 846 Dosing Interface. Touch-sensitive, high-resolution color display, simple and intuitive operation, thanks to Favorites for direct method access. With integrated Ethernet interface for direct connection to the Internet and USB interface for connecting USB printers or a USB memory stick.

Dialog languages: German, English, Chinese, French, Spanish, Portuguese, Russian, Korean, Polish and Italian.

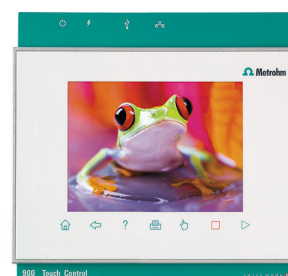

Intelligent, combined pH electrode with integrated memory chip for storing sensor data and Pt1000 temperature sensor. This electrode is particularly suitable:

- for pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- for long-term measurements




The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

Alternatively: reference electrolyte for measurements at $T > 80^\circ\text{C}$: Idrolyte, storage in Idrolyte.

iTrodes can be connected to Titrando, Ti-Touch or 913/914 meters.

To attach an 801, 804, 803 Stirrer to a Titrando, Titrino plus, Dosimat plus or 856 and 867.

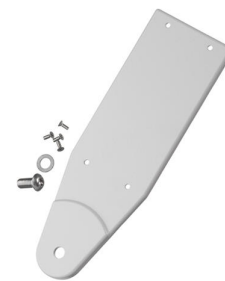

For support rods with a diameter of 10 mm.





Electrode holder for 4 electrodes and 2 buret tips



Mixed buffer solutions pH 4.00/7.00/9.00 (25 °C) in single use sachets, colourless, box of 3 x 10 x 30 mL



Maintenance kit for pH electrodes

The kit contains:

- 50 mL cleaning solution
- 50 mL 3M KCl solution
- 50 mL storage solution
- 2 Storage vessels
- Instructions for use



1 PCS

6.2621.070

5 mm hex key for IC Sample Processors



1 PCS

6.2621.130

Hexagon key 2 mm

2 mm.



1 PCS

6.6064.010

USB Memory Stick for 900 Touch Control



Optionales Zubehör

Order no.	Beschreibung	
2.801.0010	801 Stirrer Magnetic stirrer without stand for supplementing the Titrino plus, Dosimat plus, Titrandos, Sample Processors, 805 Dosimats and 780/781 pH Meters. With permanently attached cable for MSB (Metrohm Serial Bus).	
2.802.0040	802 Stirrer for 804 Ti Stand Rod stirrer including 6.1909.010 stirring propeller.	
2.804.0010	804 Ti Stand without stand rod Titration stand and controller for 802 Stirrer. The 804 Ti Stand together with the optional 802 Stirrer provides an alternative to the magnetic stirrer. Without stand.	

USB Sample Processor with one workstation and one built-in membrane pump for the automatic processing of routine samples in series with small to medium quantities. In addition to the built-in pump, an additional one (membrane or peristaltic) and up to three dosing devices for Liquid Handling tasks can be connected.

Because of the multitude of application variants, rack, stirrer, titration head, Swing Head and sample vessels must be tailored to the application and ordered separately.

The control is "stand alone" using Touch Control. The following software products can be selected for the PC control: tiamo™ titration software, MagIC Net chromatography software, viva voltammetry software, or OMNIS.



USB Sample Processor with one workstation and two built-in membrane pumps for the automatic processing of routine samples in series with small to medium quantities. In addition to the built-in pumps, up to three dosing devices for Liquid Handling tasks can be connected.

Because of the multitude of application variants, rack, stirrer, titration head, Swing Head and sample vessels must be tailored to the application and ordered separately.

The control is "stand alone" using Touch Control. The following software products can be selected for the PC control: tiamo™ titration software, MagIC Net chromatography software, viva voltammetry software, or OMNIS.



USB Sample Processor with one workstation for the automatic processing of routine samples in series with small to medium quantities. Up to two pumps (membrane or peristaltic) and three dosing devices for Liquid Handling tasks can be connected.

Because of the multitude of application variants, rack, stirrer, titration head, Swing Head and sample vessels must be tailored to the application and ordered separately.

The control is "stand alone" using Touch Control. The following software products can be selected for the PC control: tiamo™ titration software, MagIC Net chromatography software, viva voltammetry software, or OMNIS.



Stacking frame for 846 Dosing Interface, 856 Conductivity Module, 867 pH Module

For fixing the Reagent Organizer on top of the Dosing Interface



Adapter red 2 mm plug / 4 mm socket

For connecting plug B (4 mm) to 2 mm socket.



6.2103.140 Adapter black 2 mm plug / B socket 4 mm
For connecting plug B (4 mm) to 2 mm socket.



6.2151.000 Cable USB A – mini-DIN 8-pin
Controller cable

