



Getrennte Ag-Ringelektrode

6.00350.100

Diese Silberelektrode muss zusammen mit einer Referenzelektrode ($c(\text{KNO}_3) = 1 \text{ mol/L}$ als Referenzelektrolyt) verwendet werden und eignet sich für Fällungstitrationsen (Titrimittel Silbernitrat), wie z.B. von:

- Chlorid, Bromid, Iodid
- Sulfiden
- Schwefelwasserstoff
- Mercaptane
- Cyanide

Der fest eingeschmolzene Silberring ist robust gegenüber höher konzentrierten Säuren und Salzlösungen.

Lieferumfang 6.00350.100

Qt.	Order no.	Beschreibung
-----	-----------	--------------

1 PCS

6.1236.050

Sleeve with SGJ 14/12 mm

Sleeve with SGJ 14/12 mm, polyethylene.



1 PCS

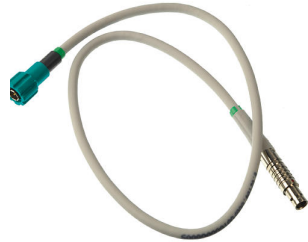
6.2008.040

Storage vessel

Together with 6.2008.050 storage vessel holder. it provides a support for the electrode on 807 Dosing Units.



Optionales Zubehör

Order no.	Beschreibung
6.02104.000	Electrode cable plug-in head G / plug P, 0.55 m Cable for connecting an analog electrode (plug-in head G) to an analog measuring module from OMNIS. 
6.0726.100	Ag/AgCl reference electrode (length 12.5 cm) Silver/silver chloride reference electrode with Double-Junction System, installation length 10 cm. The standard ground-joint 14/15 enables easy assembly, and the flexible ground-joint diaphragm, which is insensitive to contamination, can be replaced at any time. The reference electrolyte and the bridge electrolyte can be selected freely according to use, and are easy to replace. This sensor is supplied without electrolyte filling. 
6.0750.100	LL ISE reference electrode Silver / silver chloride reference electrode with double junction system. This reference electrode is well suited for: • automated applications • ion measurements • surfactant titrations The ground-joint diaphragm, which is insensitive to contamination, offers a constant and reproducible electrolyte outflow. Additionally, the reference electrolyte is gelified for even better signal stability. The sensor is delivered with $c(\text{KCl}) = 3 \text{ mol/L}$ as bridge electrolyte, which can be freely selected and exchanged as needed. 

6.2104.020

Electrode cable / 1 m / F

For connecting electrodes with Metrohm plug-in head G to Metrohm instruments (socket F).



6.2310.010

Electrolyte KNO₃-1 mol/L 250 mL

Electrolyte solution KNO₃ 1 M (reference electrolyte for combined silver electrodes)

