



Ag-Titrode

6.00430.100

Kombinierte Silberringelektrode mit einer pH-Glasmembran als Referenzelektrode.

Diese wartungsfreie Elektrode eignet sich für Fällungstitrationsen bei konstant bleibendem pH-Wert (Titriermittel Silbernitrat) von z.B.:

- Chlorid, Bromid, Iodid
- Sulfiden
- Schwefelwasserstoff
- Mercaptane
- Cyanide

Diese Elektrode wird in dest. Wasser aufbewahrt.

Je nach Anwendung empfiehlt sich die Verwendung einer Ag-Titrode mit Überzug (AgBr-, AgCl- oder Ag₂S-Überzug), die entsprechend bestellt werden kann.

Lieferumfang 6.00430.100

Qt.	Order no.	Beschreibung
-----	-----------	--------------

1 PCS

6.1236.050

Sleeve with SGJ 14/12 mm

Sleeve with SGJ 14/12 mm, polyethylene.



1 PCS

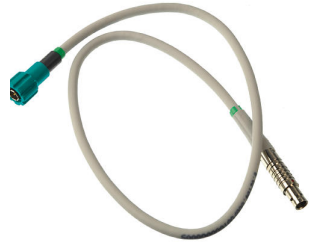
6.2008.040

Storage vessel

Together with 6.2008.050 storage vessel holder. it provides a support for the electrode on 807 Dosing Units.



Optionales Zubehör

Order no.	Beschreibung	
6.02104.000	Electrode cable plug-in head G / plug P, 0.55 m Cable for connecting an analog electrode (plug-in head G) to an analog measuring module from OMNIS.	
6.02104.010	Electrode cable plug-in head G / plug P, 1.5 m Cable for connecting an analog electrode (plug-in head G) to an analog measuring module from OMNIS.	
6.2104.020	Electrode cable / 1 m / F For connecting electrodes with Metrohm plug-in head G to Metrohm instruments (socket F).	