



Peristaltic (4-channel) pump module

2.1016.0110

Module for installation in the module plate of the Pick&Place of the OMNIS Sample Robot. This workstation is equipped with two rinsing pumps and two aspiration pumps each. They are used to clean the sensors in two Pick&Place modules and to empty the sample beaker after the analysis before it is placed back in the rack once again.

Below, the accessories are grouped into Scope of delivery and Optional accessories.
Please keep this printout at hand for ordering replacement material.
These lists may be subject to change.

Scope of delivery 2.1016.0110

Qt.	Order no.	Description
-----	-----------	-------------

2 PCS

6.01543.100

Aspiration tip automation

Aspiration tip for automated OMNIS titrations.



1 PCS

6.01803.000

Silicone tubing 8/12 2.5 m

Silicone tubing for connecting to the drainage of the Pick&Place module or pump module in order to channel any solution overflows that may occur safely into the waste canister.



1 PCS

6.01804.500

Tubing adapter for drainage tubing

Tubing adapter for drainage tubing on Pick&Place modules and pump modules in the OMNIS Sample Robot Pick&Place.




6 PCS

6.01805.000

FEP tubing connection / M6 / 66.5 cm

With light protection and protection against kinks.

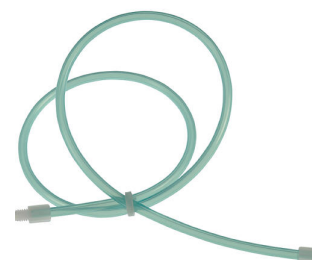

2 PCS

6.01805.010

PTFE tubing connection / M8 / 66.5 cm

With light protection and protection against kinks.



2 PCS

6.05001.010

"Peristaltic" (2-channel) pump module accessory set

The accessory set contains all of the components required for the tube connections between a rinsing pump and an aspiration pump at the "Peristaltic" pump module. In addition to preassembled tubing with Luer connectors, each of these also includes two press clamps and pump tubing.



4 PCS

6.1621.000

Container 10 L

As rinsing or waste container in automated systems.



4 PCS

6.1621.100

Cap for 10-L container

Cap for 10-L container. For easy and safe transportation of filled containers.



4 PCS

6.1828.000

PVDF connection nipple

For 6.1621.000 canister.



2 PCS

6.2709.070

Guiding sleeve

For fixation of buret and dosing tips in titration heads and electrode holders



6 PCS


6.2740.020

Spray nozzle

For the fine-spraying of the rinsing solution.



Optional accessories

Order no.	Description
6.05700.010	<div>Consumable Kit OMNIS Pump Module (2-channel)</div> <div>Two pump tubings each and associated press clamps for replacement in the OMNIS pump module. One Consumable Kit can be used for re-equipping two peristaltic pumps of the pump module.</div> <div></div>