



"Peristaltic" (2-channel) pump module accessory set

6.05001.010

The accessory set contains all of the components required for the tube connections between a rinsing pump and an aspiration pump at the "Peristaltic" pump module. In addition to preassembled tubing with Luer connectors, each of these also includes two press clamps and pump tubing.

Below, the accessories are grouped into Scope of delivery and Optional accessories.
Please keep this printout at hand for ordering replacement material.
These lists may be subject to change.

Scope of delivery 6.05001.010

Qt.	Order no.	Description
-----	-----------	-------------

2 PCS

6.01801.020

Pharm-A-line pump tubing to the peristaltic pump module

The pump tubing comes pre-equipped with adapters that permit direct installation in the peristaltic pump of the pump module. The Pharm-A-line pump tubing is used for the transport of both aqueous and nonaqueous solutions and must be operated in connection with the white press clamp in the pump module.



1 PCS

6.01805.100

PTFE tubing with Luer 60 cm

Tubing for the connection between peristaltic pump 1 or 3 in the pump module with the distributor (rinsing) on the Pick&Place module.




1 PCS

6.01805.110

PTFE tubing with Luer 80 cm

Tubing for the connection between peristaltic pump 2 or 4 in the pump module with the distributor (aspiration) on the Pick&Place module.

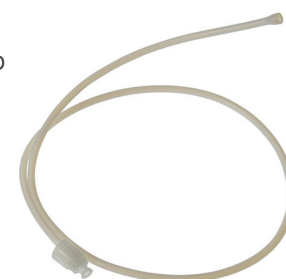

2 PCS

6.01805.120

PTFE tubing with Luer 2.5 m

Tubing for the connection between the canister and the peristaltic pump in the pump module.



2 PCS

6.02703.000

Press clamp, white, for "Peristaltic" pump module

The press clamp is used to fix and secure the pump tubing in the "Peristaltic" pump module.



1 PCS

6.02704.000

Tubing clip to workstation (10 pcs)

For fixing in place the PVC or PTFE tubing that leads to the workstation.



Optional accessories

Order no.	Description
6.1828.000	PVDF connection nipple For 6.1621.000 canister.