



IC equipment: MiPuT

6.5330.170

Accessory set for assembling a Dosino for sample pick-up (Metrohm intelligent Pick-up Injection Technique).

Below, the accessories are grouped into Scope of delivery and Optional accessories.
Please keep this printout at hand for ordering replacement material.
These lists may be subject to change.

Scope of delivery 6.5330.170

Qt.	Order no.	Description
1 PCS	6.1014.200	Metrosep I Trap 1 - 100/4.0
<p>The Metrosep I Trap 1 - 100/4.0 column is used to eliminate ionic, i.e. cationic and anionic, contaminants from aqueous solutions. Its use is especially recommended for the purification of the transfer water in combination with "MISP" or "Metrohm Inline Sample Preparation". Using the Metrosep I Trap 1 - 100/4.0 column can significantly reduce the influence of the transfer water on the system blank.</p>		



1 PCS

6.1580.120

807 Dosing Unit 2 mL without accessories

807 Dosing Unit with integrated data chip with 2 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread.



1 PCS

6.1602.160

Bottle attachment GL 45 for eluents and reagents

For eluent, reagent, and rinsing bottles, with connectors for the adsorber tube and the aspiration tubing.



1 PCS

6.1608.070

Eluent bottle / 2 L / GL 45




1 PCS

6.1618.020

Thread adapter S 40 to GL 45

Adapter for bottles with S40 thread (Merck...) to GL45.

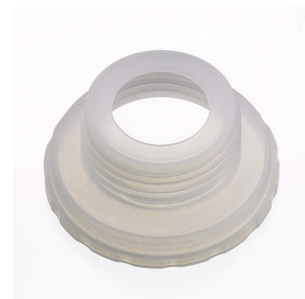

1 PCS

6.1619.000

Adsorber tube for 807 Dosing Unit



1 PCS

6.1624.000

Adapter SGJ 14 for 6.1619.XXX Adsorber tube

Adapter SGJ 14. Part of 6.1619.010 for Dilution Sample Processor.



1 PCS

6.1805.120

FEP tubing / M6 / 100 cm

With light and kink protection



1 PCS

6.1808.090

Thread adapter M8 outer / M6 inner

Outer thread M8, inner thread M6, e.g. for mounting 6.2418.XXX Air tubes in 6.1429.XXX Reaction vessels of stability measuring instruments.



1 PCS

6.1819.110

FEP aspiration tubing to canister

Aspiration tubing to canister attachment with M6 connector. For the production of eluents.



1 PCS

6.1825.290

PEEK sample loop 250 μ L

Sample loop for «Metrohm intelligent Partial loop Technique» (MiPT) and for PILS (particle into liquid sampler). With 2 Peek pressure screws.



1 PCS

6.1831.100

PEEK capillary 0.25 mm i.d., 1 m



1 PCS

6.1841.000

PEEK transfer capillary 2 mL, 5 m




1 PCS

6.2057.210

Holder for Dosino to IC instruments

Mounting part for Dosinos at the IC instrument housing.

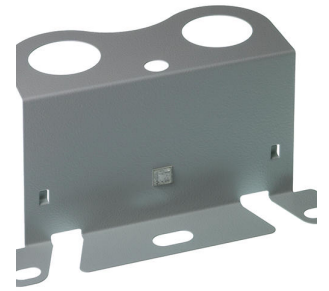

1 PCS

6.2624.100

Pick-up needle MF

Coated steel needle for metalfree pick-up sampling as well as for Eppendorf sample vials.



1 PCS

6.2744.010

Pressure screw 5x

With UNF 10/32 connection. For the connection of PEEK capillaries



1 PCS

6.2744.080

M6 thread / UNF 10/32 adapter

With M6 outer thread. For IC instruments.



1 PCS

6.2833.010

Needle/tube holder 1/16 in.

For IC Sample Processors. For 1/16 in. tubes



1 PCS

6.2841.120

Liquid Handling Station, left

Combined rinsing and dilution station for sample changers with a Swing Head. To be assembled to the left of the tower.



1 PCS

6.9920.231

Retaining plate for sample rack 159 x 2 mL + 3 x 300 mL

Plate for deflecting the vessel cover on the Professional IC Sample Processor. For 2 mL sample vessels.

