



858 Professional Sample Processor CVS con Dosino para la transferencia de muestras

2.858.0120

El 858 Professional Sample Processor CVS con Dosino para la transferencia de muestras permite el procesamiento completamente automático de grandes series de muestras de aditivos orgánicos en baños galvánicos como abrillantadores y supresores mediante voltamperometría de redisolución cíclica (CVS, por sus siglas en inglés) junto con los aparatos Professional VA/CVS y el software viva. El sistema se controla directamente a través del software viva y permite el cambio automático de métodos entre las distintas clases de determinación. La gradilla de muestras suministrada tiene capacidad para 56 muestras de un máximo de 11 mL para la determinación de supresores y 28 muestras de un máximo de 50 mL para la determinación de abrillantadores. Se utiliza un Dosino para la transferencia de muestras en el vaso de medida. De este modo, se puede reducir la cantidad de muestra necesaria y medirla automáticamente.

El 858 Professional Sample Processor CVS se entrega con Dosino para la transferencia de muestras con el brazo giratorio montado, incluidos los accesorios, el Dosino, la gradilla de muestras y el cable de conexión para el montaje de un puesto completo de trabajo, pero sin equipo de lavado.

Partes/accesorios 2.858.0120


Qt.	Order no.	Descripción
1 PCS	1.800.0040	800 Dosino (cable length 0.65 m)
Drive with write/read hardware for intelligent dosing units. With permanently attached cable (length 0.65 m).		
		
1 PCS	1.858.0010	858 Professional Sample Processor
The 858 Professional Sample Processor processes samples from 500 µL to 500 mL. The sample transfer takes place either by means of a peristaltic pump on the 850 Professional IC system or with an 800 Dosino.		
		
1 PCS	6.1546.030	Piston tongs
For pistons for 807 Dosing Units or cylinder units OMNIS.		
		
1 PCS	6.1602.115	Bottle attachment / S 40 / green
For direct dosing from reagent bottles with S 40 thread (Merck...).		
		

1 PCS

6.1618.020

Thread adapter S 40 to GL 45

Adapter for bottles with S40 thread (Merck...) to GL45.

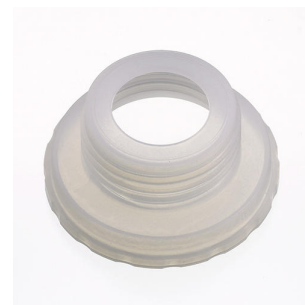

1 PCS

6.1803.020

PTFE capillary 0.97 mm i.d. / 5 m

For all IC instruments



1 PCS

6.1805.530

FEP tubing / M6 / 2 m

With light and kink protection



1 PCS

6.1808.280

Adapter 807 Dosing Unit Port 4, M6 inner

Adapter for using Port 4 of the 807 dosing Unit for connection with an FEP tubing / M6.



1 PCS

6.1815.010

Spiral band 0.5 m

For holding together different cables or tubing.



1 PCS

6.1819.010

PTFE tube

As gas inlet in washing bottle of VA Stand



1 PCS

6.1831.020

PEEK capillary 0.75 mm i.d., 22 cm

For all IC components.



1 PCS

6.1835.050

PEEK sample tube

For sample changers for transferring larger-sized sample volumes, e.g. for brightener determination.




1 PCS

6.2041.450

Sample rack 56 x 11 mL + 56 x 50 mL

Multi-row sample rack for 56 samples in 11 mL sample vessels (6.2743.050) and 56 samples in sample vessels of 50 mL (6.2747.010).

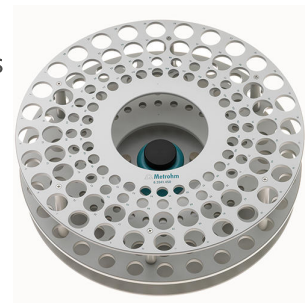

1 PCS

6.2057.040

Holder for Dosino for Sample Processor

Holder for Dosino for the 838 / 858 IC Sample Processor and the 838 Advanced VA Sample Processor



1 PCS

6.2151.000

Cable USB A – mini-DIN 8-pin

Controller cable



1 PCS

6.2621.080

Capillary cutter

Cutting tool for simple cutting of plastic capillaries, 1/16" and 1/8", including 5 replacement blades.



1 PCS

6.2621.090

1/2 in. wrench

For 1/2 in. screws. For IC Sample Processors



1 PCS

6.2621.100

3 mm hex key for IC Sample Processors



1 PCS

6.2621.130

Hexagon key 2 mm

2 mm.



1 PCS

6.2739.000

Wrench

For tightening connectors



1 PCS

6.2743.057

Sample vessel 11 mL

200 pieces. For IC Sample Processors and VA Autosamplers.

The IC sample vessels are made of polypropylene (PP) plastic for health-care applications. Production and packaging both take place under sterile and clean room conditions (Clean Room Class 7, in accordance with DIN ISO 14644). This guarantees minimal degrees of contamination. The sample vessels are subjected to random checks for leachable anions and cations.



1 PCS

6.2744.010

Pressure screw 5x

With UNF 10/32 connection. For the connection of PEEK capillaries



2 PCS

6.2744.080

M6 thread / UNF 10/32 adapter

With M6 outer thread. For IC instruments.



4 PCS

6.2747.010

Sample vessel 50 mL

With blue screw cap and mL graduation. Set of 25 pieces.




1 PCS

6.2751.110

Safety shield

Safety shield for the Advanced IC/VA Sample Processors.

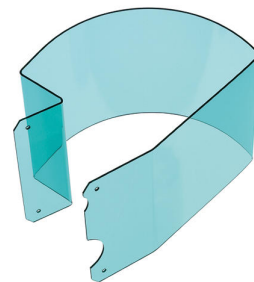

1 PCS

6.3032.150

807 Dosing Unit 5 mL

807 Dosing Unit with integrated data chip with 5 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.

