



Ag-Titrode

6.00430.100

Electrodo de anillo de plata combinado con una membrana de vidrio de pH como electrodo de referencia.

Este electrodo no precisa mantenimiento y es apto para titulaciones por precipitación con un valor de pH constante (reactivo de titulación: nitrato de plata), por ejemplo, a partir de:

- Cloruro, bromuro, yoduro
- Sulfuros
- Ácido sulfhídrico
- Mercaptano
- Cianuro

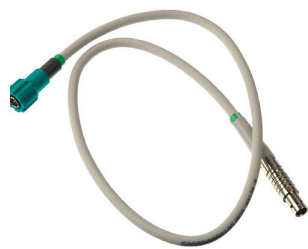
Este electrodo se conserva en agua destilada.

Dependiendo de la aplicación, se recomienda el uso de un Ag-Titrode con recubrimiento (de AgBr, AgCl o Ag₂S), que se puede incluir en el pedido.

Partes/accesorios 6.00430.100

Qt.	Order no.	Descripción
1 PCS	6.1236.050	Sleeve with SGJ 14/12 mm Sleeve with SGJ 14/12 mm, polyethylene.
		
1 PCS	6.2008.040	Storage vessel Together with 6.2008.050 storage vessel holder. it provides a support for the electrode on 807 Dosing Units.
		

Accesorios opcionales

Order no.	Descripción	
6.02104.000	Electrode cable plug-in head G / plug P, 0.55 m Cable for connecting an analog electrode (plug-in head G) to an analog measuring module from OMNIS.	
6.02104.010	Electrode cable plug-in head G / plug P, 1.5 m Cable for connecting an analog electrode (plug-in head G) to an analog measuring module from OMNIS.	
6.2104.020	Electrode cable / 1 m / F For connecting electrodes with Metrohm plug-in head G to Metrohm instruments (socket F).	