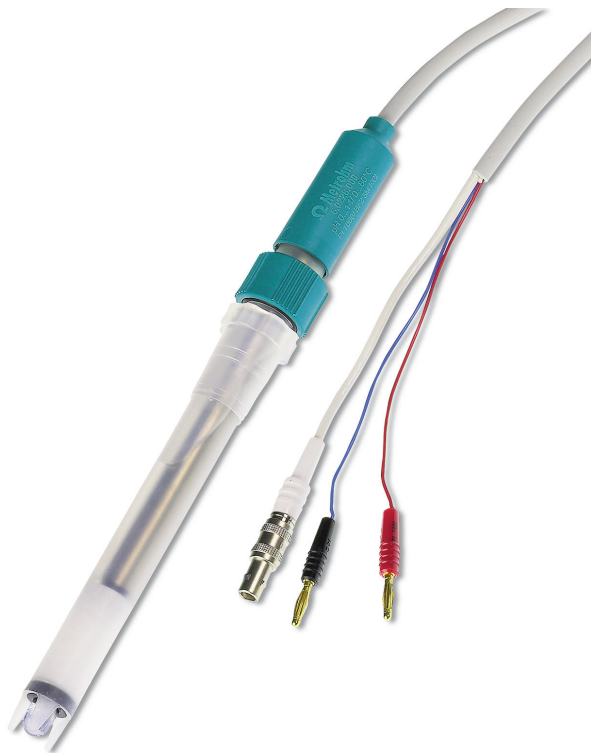




Solitrode con Pt1000 (cable fijo de 1,2 m)

6.0228.000

Electrodo pH combinado con sensor de temperatura Pt1000 integrado y cable fijo (1,2 m). Este electrodo es apto para las medidas de pH rutinarias en soluciones libres de precipitados, proteínas y sulfuros. Gracias al eje de plástico robusto/irrompible de polipropileno y la protección de la membrana de vidrio contra impactos, el electrodo es mecánicamente muy resistente.

Electrolito de referencia: $c(\text{KCl}) = 3 \text{ mol/L}$, almacenamiento en una solución de conservación.

Partes/accesorios 6.0228.000

Qt.	Order no.	Descripción
1 PCS	6.1236.050	Sleeve with SGJ 14/12 mm
		Sleeve with SGJ 14/12 mm, polyethylene.



Together with 6.2008.050 storage vessel holder, it provides a support for the electrode on 807 Dosing Units.



Accesorios opcionales

Order no.	Descripción	
6.2308.020	Electrolyte 3 mol/L KCl (250 mL) Electrolyte solution $c(\text{KCl}) = 3 \text{ mol/L}$, (for Ag/AgCl reference systems)	
6.2323.000	Storage solution Storage solution for all combined pH glass electrodes with reference electrolyte $c(\text{KCl}) = 3 \text{ mol/L}$	
6.2325.000	pHit kit Maintenance kit for pH electrodes The kit contains: • 50 mL cleaning solution • 50 mL 3M KCl solution • 50 mL storage solution • 2 Storage vessels • Instructions for use	

Reliable measuring results over long periods of time can be guaranteed only if the pH glass membrane and the diaphragm receive preventive and regular care. Cleaning by etching with toxic chemicals or applying mechanical treatment to the diaphragm is not only complicated and expensive, it also accelerates the ageing of the pH glass electrode.

The cleaning solution was developed for easy and gentle cleaning of pH glass electrodes. Regular use can considerably prolong their service life.



This cleaning solution is also part of the pHit kit.
