



Metrosep A Supp 7 - 250/4,0

6.1006.630

Se sospecha que los subproductos del tratamiento del agua (subproductos de la desinfección) no solo pueden resultar nocivos para la salud, sino también ser cancerígenos. Por eso, los oxohalogenuros han sido objeto de muchos estudios y normas (p. ej. EPA 300.1 parte B, EPA 317.0, EPA 326.0). En particular, se trata del bromato que se forma en la ozonificación de agua potable a partir de bromuro. La Metrosep A Supp 7 - 250/4,0 es una columna de separación de alto rendimiento para la determinación paralela de aniones estándar, oxohalogenuros y ácido dicloroacético. Con esta columna es posible determinar estos iones hasta la gama de $\mu\text{g/L}$ inferior de forma segura y precisa. La alta sensibilidad de detección se logra usando el polímero de alcohol polivinílico de 5 μm con el que pueden alcanzarse números de platos extremadamente elevados y con ellos excelentes propiedades de separación y detección. Además, es posible adaptar la separación a las exigencias específicas de la aplicación modificando la temperatura.

Accesorios opcionales

Order no.	Descripción
-----------	-------------

The Metrosep A Supp 5 Guard/4.0 reliably protects the Metrosep A Supp 5 and 7 IC anion columns against contamination from the sample or the eluent.

It contains the same separation material as the Metrosep A Supp 5, is also made of PEEK, and is screwed directly onto the respective separation column with virtually no dead volume ("On Column Guard System"). The guard column prolongs the service life of the analytical column, with practically no influence on its chromatographic separating efficiency. The economical price and simple handling make using the A Supp 5 Guard/4.0 highly recommended.




The Metrosep A Supp 5 S-Guard/4.0 reliably protects the Metrosep A Supp 5 and 7 IC anion columns against contaminations from the sample or eluent.



The Metrosep RP 2 Guard/3.5 is a guard column for universal use. It reliably protects the analytical separation column against contaminations. It reliably protects the analytical separation column against contamination, removing the smallest particles, traces of iron oxide and bacteria. The Metrosep RP 2 Guard/3.5 helps to reduce costs; its filter disk can be replaced very easily.




Do not use for 6.1011.030.

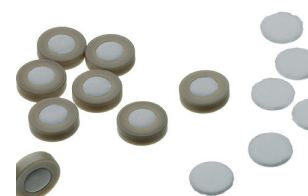

6.1011.130

Spare filter for RP 2 Guard/3.5, 10 pieces

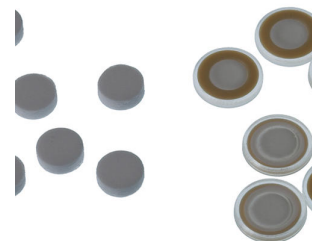

6.1031.500

Metrosep A Supp 16 Guard/4.0

The Metrosep A Supp 16 Guard/4.0 reliably protects the Metrosep A Supp 16 analytical separation columns against contamination. Thanks to the "On Column Guard System", the guard column is very easy to handle. The guard column screws easily onto the analytical column. No tools are required.



6.1031.510

Metrosep A Supp 16 S-Guard/4.0

The Metrosep A Supp 16 S-Guard/4.0 reliably protects the Metrosep A Supp 16 analytical separation columns against contamination. The guard column is connected to the separation column simply using capillary connections. No tools are required.

