



OMNIS Titrator sans agitateur, sans licence fonctionnelle

2.1001.0010

OMNIS Titrator, innovant, modulaire, potentiométrique. Grâce à la technologie 3S de l'adaptateur Liquid Adapter, la manipulation des produits chimiques est plus sûre que jamais. Avec des modules de mesure et des unités de cylindre, le titreur peut être librement configuré et il est possible au besoin d'y ajouter un agitateur. Grâce à différentes licences fonctionnelles du logiciel, des modes de mesure et des fonctionnalités variées sont possibles.

- Commande via un PC ou un réseau local
- Possibilité de connecter jusqu'à quatre autres modules de titrage ou de dosage pour d'autres applications ou solutions auxiliaires
- Possibilité d'y ajouter un agitateur magnétique et/ou un agitateur à hélice
- Différentes tailles de cylindre disponibles : 5, 10, 20 ou 50 mL
- Liquid Adapter avec la technologie 3S : manipulation de produits chimiques plus sûre et transfert automatique des données des réactifs

Modes de mesure et options logicielles :

- Titrage à point final : licence fonctionnelle « Basic »

- Titration à point final et à point d'équivalence (monotone/dynamique) : licence fonctionnelle « Advanced »
- Titration à point final et à point d'équivalence (monotone/dynamique) avec titration en parallèle : licence fonctionnelle « Professional »
- Licences optionnelles pour d'autres modes de mesure, de titration et d'exploitation (par exemple titration volumétrique Karl-Fischer, titration conductimétrique ou thermométrique, coulométrie, détermination de l'indice de brome)

Contenu de la livraison 2.1001.0010

Qt.	Order no.	Description	
1 PCS	6.01600.010	OMNIS Liquid Adapter	
		Adapter for connecting reagent bottles with 3S technology to an OMNIS Titrator, Titration Module or Dosing Module.	
1 PCS	6.01601.000	Bottle cap, multi-use	
		Reusable bottle cap with 3S technology for connecting reagent bottles to an OMNIS Titrator, Titration Module or Dosing Module.	
1 PCS	6.02004.010	Bottle holder	
		Bottle holder for OMNIS Titrator, Titration Module and Dosing Module.	

1 PCS

6.02005.010

Electrode Parking Station

Electrode Parking Station for electrode storage on an OMNIS titration system.



1 PCS

6.1608.020

Glass bottle / 1000 mL / GL 45

Bottle for auxiliary solutions.



Accessoires optionnels

Order no.	Description	
2.1002.0010	OMNIS Titration Module without stirrer Titration module for connection to an OMNIS Titrator for extending the system to include an additional buret for titration/dosing and up to two measuring modules (analog or digital). Can be supplemented with one magnetic stirrer or rod stirrer for use as separate titration stand. Freely selectable cylinder unit with 5, 10, 20 or 50 mL.	
2.1002.0110	OMNIS Titration Module with magnetic stirrer Titration module with magnetic stirrer for connection to an OMNIS Titrator for extending the system to include an additional titration stand, a buret for titration/dosing and up to two measuring modules (analog or digital). Freely selectable cylinder unit with 5, 10, 20 or 50 mL.	
2.1003.0010	OMNIS Dosing Module without stirrer Dosing module for connection to an OMNIS Titrator for extending the system to include an additional buret for titration/dosing. Can be supplemented with one magnetic stirrer or rod stirrer for use as separate titration stand. Freely selectable cylinder unit with 5, 10, 20 or 50 mL.	
2.1003.0110	OMNIS Dosing Module with magnetic stirrer Dosing module with magnetic stirrer for connection to an OMNIS Titrator for extending the system to include an additional titration stand and a buret for titration/dosing. Freely selectable cylinder unit with 5, 10, 20 or 50 mL.	

Rod stirrer for connection to one OMNIS Titrator (cable length 1.1 m).



6.00200.300 dUnitrode with Pt1000

Digital, combined pH electrode for OMNIS with integrated Pt1000 temperature sensor. This electrode is particularly suitable for :

- pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- long-term measurements

The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

Alternatively: Reference electrolyte for measurements at $T > 80^\circ\text{C}$: Idrolyte, storage in Idrolyte.

dTrodes can be used on OMNIS Titrators.



6.00201.300 dEcotrode Plus

Digital, combined pH electrode for OMNIS.

The electrode is suitable for aqueous acid/base titrations.

The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

dTrodes can be used on OMNIS Titrators.



Digital, combined pH electrode for OMNIS, with integrated Pt1000 temperature sensor, for pH measurements/titrations in ion-deficient aqueous media (e.g., drinking water, process water). This electrode shows a very short response time in these samples.

The fixed ground-joint diaphragm is insensitive to contamination.

When $c(\text{KCl}) = 3 \text{ mol/L}$ is used as bridge electrolyte, storage in storage solution is recommended.

The bridge electrolyte can be easily replaced with a chloride-free electrolyte (e.g., potassium nitrate $c(\text{KNO}_3) = 1 \text{ mol/L}$ (6.2310.010)).

Storage in the bridge electrolyte used.

dTodes can be used on OMNIS Titrators.



Digital, combined pH electrode for OMNIS for all non-aqueous acid/base titrations. The glass membrane is optimized for poorly conducting solutions and thanks to the flexible ground-joint diaphragm, the electrode is suitable for contaminated samples.

This electrode can be used with non-aqueous reference electrolytes (lithium chloride or tetraethylammonium bromide).

Storage in corresponding reference electrolyte.

dTodes can be used on OMNIS Titrators.



Digital, combined pH electrode for OMNIS with double-junction system. This electrode is suitable for pH measurements/titrations of samples

- that contaminate the reference system of the electrode (e.g., galvanic baths, samples containing sulfides)
- for which potassium chloride $c(\text{KCl}) = 3 \text{ mol/L}$ cannot be used as reference electrolyte (e.g., reaction of potassium or chloride with sample)

The electrode is equipped with a flexible ground-joint diaphragm which is insensitive to contamination and which can be exchanged when needed.

When $c(\text{KCl}) = 3 \text{ mol/L}$ is used as bridge electrolyte, storage in storage solution is recommended.

The bridge electrolyte can be easily replaced with a suitable electrolyte (e.g., potassium nitrate $c(\text{KNO}_3) = 1 \text{ mol/L}$ (6.2310.010)). Storage in the bridge electrolyte used.

dTodes can be used on OMNIS Titrators.



Digital, combined platinum ring electrode for OMNIS with a pH glass membrane as reference electrode.

This maintenance-free electrode is suitable for redox titrations when the pH value remains constant, e.g.:

- Iodometry
- Chromatometry
- Cerimetry
- Permanganometry

This electrode is stored in distilled water.

dTodes can be used on OMNIS Titrators.



Digital, combined silver ring electrode for OMNIS with fixed ground-joint diaphragm.

This electrode is suitable for precipitation titrations where the pH value varies (titrant silver nitrate) of, e.g.

- Chloride, Bromide, Iodide
- Sulfides
- Hydrogen sulfide
- Mercaptans
- Cyanides

The fixed ground-joint diaphragm is insensitive to contamination. $c(\text{KNO}_3) = 1 \text{ mol/L}$ is used as reference electrolyte and for storage. dTrodes can be used on OMNIS Titrators.



Digital, combined platinum ring electrode for OMNIS with fixed ground-joint diaphragm.

This electrode is suitable for redox titrations when the pH value varies, e.g.:

- Oxygen content according to Winkler
- Determination of hydrogen peroxide with KMnO_4
- Diazotization titrations

The fixed ground-joint diaphragm is insensitive to contamination. $c(\text{KCl}) = 3 \text{ mol/L}$ is used as reference electrolyte and for storage. dTrodes can be used on OMNIS Titrators.



Digital, combined calcium-selective electrode for OMNIS.

This ISE is suitable for:

- Ion measurements of Ca^{2+} ($1 \cdot 10^{-7}$ to 1 mol/L) in aqueous solutions
- Complexometric (back) titrations (e.g., determination of water hardness)



Thanks to a robust/break-proof plastic shaft made of propylene and an impact protection for the polymer membrane, this sensor is mechanically very resistant.

The reference electrolyte used is $c(\text{KCl}) = 3 \text{ mol/L}$.

dTrodes can be used on OMNIS Titrators.

6.01503.120 OMNIS 2 mL cylinder unit without accessories

Intelligent 2 mL cylinder unit for an OMNIS Titrator, Titration Module or Dosing Module. Without dosing tubing and antidiffusion tip.



6.01503.150 OMNIS 5 mL cylinder unit without accessories

Intelligent 5 mL cylinder unit for an OMNIS Titrator, Titration Module or Dosing Module. Without dosing tubing and antidiffusion tip.



6.01503.210 OMNIS 10 mL cylinder unit without accessories

Intelligent 10 mL cylinder unit for an OMNIS Titrator, Titration Module or Dosing Module. Without dosing tubing and antidiffusion tip.



6.01503.220 OMNIS 20 mL cylinder unit without accessories

Intelligent 20 mL cylinder unit for an OMNIS Titrator, Titration Module or Dosing Module. Without dosing tubing and antidiffusion tip.



6.01503.250 OMNIS 50 mL cylinder unit without accessories

Intelligent 50 mL cylinder unit for an OMNIS Titrator, Titration Module or Dosing Module. Without dosing tubing and antidiffusion tip.



6.01900.010 Stirring propeller 30 mm ETFE

Stirring propeller 30 mm ETFE with adjustable height for use with Rod Stirrer "Titrator" or Rod Stirrer "Sample Robot"



6.01900.030 Stirring propeller 20 mm ETFE

Stirring propeller, 20 mm, with adjustable height for use with Rod Stirrer "Titrator" or Rod Stirrer "Sample Robot"



6.02100.010 Measuring module digital

Digital measurement channel for one OMNIS Titrator or Titration Module for the connection of digital "dTrodes" electrodes.



6.02101.011 Measuring module analog

Analog measurement channel for one OMNIS Titrator or Titration Module for the connection of analog electrodes. Compatible with OMNIS software version 2.18 or higher.



6.02102.010 Cable MDL PL/SO 0.5 m

Connecting cable for OMNIS instruments with MDL connector (0.5 m)



6.02102.020 Cable MDL PL/SO 1 m

Connecting cable for OMNIS instruments with MDL connector (1 m)



6.02102.030 Cable MDL PL/SO 2 m

Connecting cable for OMNIS instruments with MDL connector (2 m)



6.02102.040 Cable MDL PL/SO 5 m

Connecting cable for OMNIS instruments with MDL connector (5 m)




6.02102.050 Cable MDL PL/SO 10 m

Connecting cable for OMNIS instruments with MDL connector (10 m)



6.02104.000 Electrode cable plug-in head G / plug P, 0.55 m

Cable for connecting an analog electrode (plug-in head G) to an analog measuring module from OMNIS.

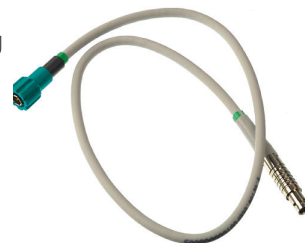

6.02104.010 Electrode cable plug-in head G / plug P, 1.5 m

Cable for connecting an analog electrode (plug-in head G) to an analog measuring module from OMNIS.




6.02104.020 Electrode cable plug-in head G (temp.) / plug P, 0.55 m

Cable for connecting a temperature sensor (head G) to an analog measuring module from OMNIS.



6.02104.030 Electrode cable plug-in head G (temp.) / plug P, 1.5 m

Cable for connecting a temperature sensor (head G) to an analog measuring module from OMNIS.

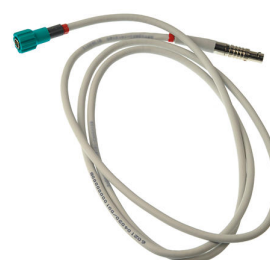

6.02104.040 Electrode cable plug-in head G (pol.) / plug P, 0.55 m

Cable for connecting an analog polarizable electrode (plug-in head G) to an analog measuring module from OMNIS.



6.02104.050 Electrode cable plug-in head G (pol.) / plug P, 1.5 m

Cable for connecting an analog polarizable electrode (plug-in head G) to an analog measuring module from OMNIS.



6.02104.300 Electrode cable plug-in head Q / plug P, 0.55 m

Cable for connecting a dTrobe (plug-in head Q) to the digital measuring module from OMNIS.



6.02104.310 Electrode cable plug-in head Q / plug P, 1.5 m

Cable for connecting a dTrobe (plug-in head Q) to the digital measuring module from OMNIS.



6.02104.600 Electrode cable plug-in head U / plug P, 0.55 m

Cable for connecting an electrode with plug-in head U to an analog measuring module from OMNIS.




6.02104.610 Electrode cable plug-in head U / plug P, 1.5 m

Cable for connecting an electrode with plug-in head U to an analog measuring module from OMNIS.




6.02106.000 Liquid Adapter extension cable, 0.55 m

Extension cable (0.55 m in length) for connecting an OMNIS Liquid Adapter to an OMNIS Titrator, Titration Module or Dosing Module.

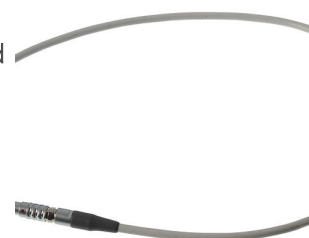

6.02106.010 Liquid Adapter extension cable, 2 m

Extension cable (2 m in length) for connecting an OMNIS Liquid Adapter to an OMNIS Titrator, Titration Module or Dosing Module.

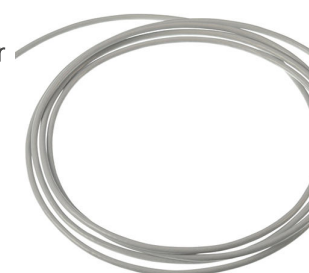

6.03001.120 Cylinder unit OMNIS 2 mL

Intelligent 2 mL cylinder unit for an OMNIS Titrator, Titration Module or Dosing Module. Includes dosing tubing and antidiffusion tip.



6.03001.150 Cylinder unit OMNIS 5 mL

Intelligent 5 mL cylinder unit for an OMNIS Titrator, Titration Module or Dosing Module. Includes dosing tubing and antidiffusion tip.



6.03001.210 Cylinder unit OMNIS 10 mL

Intelligent 10 mL cylinder unit for an OMNIS Titrator, Titration Module or Dosing Module. Includes dosing tubing and antidiffusion tip.



6.03001.220 Cylinder unit OMNIS 20 mL

Intelligent 20 mL cylinder unit for an OMNIS Titrator, Titration Module or Dosing Module. Includes dosing tubing and antidiffusion tip.



6.03001.250 Cylinder unit OMNIS 50 mL

Intelligent 50 mL cylinder unit for an OMNIS Titrator, Titration Module or Dosing Module. Includes dosing tubing and antidiffusion tip.



6.05609.000 KFT equipment for OMNIS Titrator

Complete accessory set for volumetric Karl Fischer titration.



6.05609.010 KFT equipment for OMNIS Titrator with Polytron

Complete accessory set for volumetric Karl Fischer titration with Polytron.



6.06001.100 Function license DET

"Dynamic equivalence point titration" function license for one OMNIS titration system.



6.06001.110 Function license MET

"Monotonic equivalence point titration" function license for one OMNIS titration system.



6.06001.160 Function license KFT with conditioning

Function license for volumetric Karl Fischer titration with conditioning for an OMNIS titration system.



6.06001.320 Function license for 5 parallel titrations

Function license for a maximum of five parallel titrations using one OMNIS Titrator with additional Titration Modules or Dosing Modules.



6.06001.900 Function license Basic

Function license "Basic" for one OMNIS Titrator



- Liquid Handling (LQH)
- MEAS U/T/pH
- MEAS CONC (direct measurement and standard addition)
- Calibration (Cal)
- Endpoint titration without conditioning (SET)
- Titration with OMNIS Titration Module or OMNIS Dosing Module
- Parallel dosing

Function license "Advanced" for one OMNIS Titrator



- Liquid Handling (LQH)
- MEAS U/T/pH
- MEAS CONC (direct measurement and standard addition)
- Calibration (Cal)
- Endpoint titration without conditioning (SET)
- Monotonic equivalence point titration (MET)
- Dynamic equivalence point titration (DET)
- Titration with OMNIS Titration Module or OMNIS Dosing Module
- Parallel dosing

Function license "Professional" for one OMNIS Titrator



- Liquid Handling (LQH)
- MEAS U/T/pH
- MEAS CONC (direct measurement and standard addition)
- Calibration (Cal)
- Endpoint titration without conditioning (SET)
- Monotonic equivalence point titration (MET)
- Dynamic equivalence point titration (DET)
- Titration with OMNIS Titration Module or OMNIS Dosing Module
- Parallel dosing
- Parallel titration