



913 pH Meter

2.913.0010

Appareil de mesure du pH bicanal portable pour mesurer le pH / la tension en mV et la température. Avec cet appareil de mesure fonctionnant sur batterie, vous êtes parfaitement équipé pour des mesures sur le terrain.

- pH-mètre portable avec jeu de batteries incorporé et deux entrées de mesure du pH isolées électriquement
- Entrée de mesure du pH analogique pour des électrodes pH standard de Metrohm
- Entrée de mesure du pH numérique pour des électrodes pH intelligentes de Metrohm
- Boîtier robuste, imperméable à l'eau et protégé contre la poussière (IP67), idéal pour une utilisation dans des conditions difficiles en extérieur et en laboratoire
- Écran d'affichage couleur à LCD rétro-éclairé pour une lecture facile des résultats
- Interface USB pour une exportation aisée des données sur ordinateur ou imprimante
- Mémoire interne de grande capacité (10 000 blocs de données)
- Modes utilisateur et expert protégés par PIN afin d'éviter toute modification accidentelle des paramètres
- Impression et exportation des données conformes aux BPL avec ID utilisateur et horodatage

Contenu de la livraison 2.913.0010

Qt.	Order no.	Description
1 PCS	1.913.0010	913 pH Meter
pH/mV and temperature measurement for routine use in the laboratory and on the road. pH or mV and temperature can be measured in parallel with one analog and one intelligent sensor and output to a large colour display. Important information such as charge state, user, IDs can be clearly seen at a glance. A PIN-protected expert mode protects against unintentional changes of different parameters. The meter is furnished with an accumulator for mobile use that can be charged practically anywhere. It naturally also satisfies the requirements of IP67. A stand plate allows the mobile meter to be easily converted into a laboratory meter and vice versa. Very large measured value memory (10,000 data sets) and USB interface (GLP-compliant printout or data export with optional management of the data in tiBase) offer professional data handling.		
		
1 PCS	6.2008.060	Holder for electrode storage vessel
Practical holder for fastening the electrode storage vessels to the 912, 913 or 914 pH/Conductometers.		
		
1 PCS	6.2050.010	Carrying strap for 912/913/914
Carrying strap for 912/913/914 meters		
		

1 PCS

6.2151.100

Adapter USB MINI (OTG) - USB A

For connecting USB instruments.



1 PCS

6.2151.110

Metrohm USB Mini B cable (OTG) - USB A, 1.8 m

For connecting USB instruments.



1 PCS

6.2166.100

USB power supply unit 5.25 V / 1.53 A

USB power supply unit for 912 / 913 / 914

Efficiency Level VI



Accessoires optionnels

Order no.	Description
2.142.0100	Custom Q3X thermal printer Compact printer with USB interface for • 900 Touch Control • 915 KF Ti-Touch • 916 Ti-Touch • 917 Coulometer • 877 / 848 Titrino plus • 865 / 876 Dosimat plus • 91X Meter (cable 6.2151.140) • Eco Dosimat / Titrator • 862 Compact Titrosampler • 870 KF Titrino plus • 899 Coulometer Paper width 60 mm (40 characters). Including 6.2151.120 USB cable.
2.854.0010	854 iConnect The 854 iConnect – Electrode cable and measuring amplifier for intelligent "iTrodes" electrodes (cable length 150 cm).
6.00226.600	Spearhead electrode with Pt1000 Maintenance-free combined pH electrode (gel electrolyte) for piercing measurements of all types (e.g., with cheese, meat, dough) with integrated Pt1000 temperature sensor. The electrode is stored in saturated potassium chloride solution $c(\text{KCl}) = \text{sat.}$ (6.2308.000) and is not suitable for low-ion solutions. The ageing indicator gives early indication of when the electrode needs to be replaced.



Combined pH electrode for measurements in very small sample volumes (>50 µL) and biological samples.

Idrolyte (6.2308.040) is used as reference electrolyte and storage solution.



Combined pH electrode with integrated Pt1000 temperature sensor and fixed cable (1.2 m). This electrode is suitable for routine pH measurements in solutions that do not contain precipitates, proteins, or sulfides. This electrode is mechanically resistant thanks to the robust /unbreakable plastic shaft made of polypropylene and impact protection for the glass membrane.

Reference electrolyte: c(KCl) = 3 mol/L, storage in storage solution.



Combined pH electrode with integrated temperature sensor (Pt1000) and fixed cable conforming to IP67 (1.2 m).

This electrode is suitable for an introduction to GLP-compliant pH measurements in solutions that do not contain precipitates, proteins, or sulfides. This electrode is mechanically resistant thanks to the robust /unbreakable plastic shaft made of polypropylene and impact protection for the glass membrane.

Additionally, this electrode offers a waterproof plug in accordance with IP67 for mobile use with Metrohm pH meters.

Reference electrolyte: c(KCl) = 3 mol/L, storage in storage solution.




Combined pH electrode with integrated Pt1000 temperature sensor and fixed cable (1.2 m, diameter banana plug 2 mm). This electrode is particularly suitable

- for pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- for long-term measurements

The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

Alternatively: reference electrolyte for measurements at $T > 80^\circ\text{C}$: Idrolyte, storage in Idrolyte.

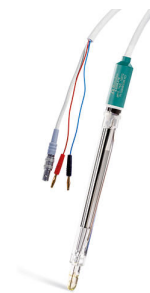

Intelligent, combined pH electrode with integrated memory chip for storing sensor data and Pt1000 temperature sensor. This electrode is particularly suitable:

- for pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- for long-term measurements

The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

Alternatively: reference electrolyte for measurements at $T > 80^\circ\text{C}$: Idrolyte, storage in Idrolyte.

iTrodes can be connected to Titrand, Ti-Touch or 913/914 meters.



Stand plate for converting a mobile 912/913/914 pH/Conductometer into a laboratory meter.



6.2104.600

Electrode cable for plug in head U/plug F, 2x2 mm B, 1m

For connecting electrodes with Metrohm plug-in head U to Metrohm instruments (socket F).



6.2151.140

Y cable USB A St - USB B St - Mini B St

Y cable for connection of a USB printer to the pH/Conductometers 912 / 913 / 914. This cable allows a printer and the power supply unit to be connected to the measuring instrument at the same time.



6.2166.500

12 V USB adapter for 912 / 913 / 914 pH/Conductometer

12 V USB adapter for 912 / 913 / 914 pH/Conductometer.



6.2307.230

Buffer solutions pH 4, 7 and 9

Mixed buffer solutions pH 4.00/7.00/9.00 (25 °C) in single use sachets, colourless, box of 3 x 10 x 30 mL



6.2313.000

Electrolyte 3 mol/L KCl (1000 mL)

Electrolyte solution $c(\text{KCl}) = 3 \text{ mol/L}$ (for Ag/AgCl reference systems)



6.2325.000

pHit kit

Maintenance kit for pH electrodes

The kit contains:

- 50 mL cleaning solution
- 50 mL 3M KCl solution
- 50 mL storage solution
- 2 Storage vessels
- Instructions for use

