



Metrosep A Supp 7 - 250/4.0

6.1006.630

Les produits secondaires du traitement de l'eau (disinfection by-products) sont soupçonnés d'être non seulement nocifs pour la santé, mais aussi cancérigènes. C'est la raison pour laquelle les oxyhalogénures sont devenus l'objet de nombreuses recherches et de normes (p. ex. EPA 300.1 partie B, EPA 317.0, EPA 326.0). Il s'agit en premier lieu du bromate, produit du bromure lors de l'ozonisation de l'eau potable. La colonne Metrosep A Supp 7 - 250/4.0 est une colonne de séparation haute performance pour la détermination en parallèle des anions standard, des oxyhalogénures et de l'acide dichloracétique. Avec cette colonne, ces ions peuvent être déterminés jusque dans la gamme inférieure des µg/L de manière sûre et précise. La sensibilité de détection élevée est atteinte par l'utilisation d'un polymère à base d'alcool polyvinylique de 5 µm, permettant d'obtenir un nombre de plateaux très élevé et, par conséquent, d'excellentes propriétés de séparation et de détection. En plus, par modification de la température, il est possible d'adapter la séparation aux exigences spécifiques de l'application.

Accessoires optionnels

Order no.	Description
-----------	-------------

The Metrosep A Supp 5 Guard/4.0 reliably protects the Metrosep A Supp 5 and 7 IC anion columns against contamination from the sample or the eluent.

It contains the same separation material as the Metrosep A Supp 5, is also made of PEEK, and is screwed directly onto the respective separation column with virtually no dead volume ("On Column Guard System"). The guard column prolongs the service life of the analytical column, with practically no influence on its chromatographic separating efficiency. The economical price and simple handling make using the A Supp 5 Guard/4.0 highly recommended.



The Metrosep A Supp 5 S-Guard/4.0 reliably protects the Metrosep A Supp 5 and 7 IC anion columns against contaminations from the sample or eluent.



The Metrosep RP 2 Guard/3.5 is a guard column for universal use. It reliably protects the analytical separation column against contaminations. It reliably protects the analytical separation column against contamination, removing the smallest particles, traces of iron oxide and bacteria. The Metrosep RP 2 Guard/3.5 helps to reduce costs; its filter disk can be replaced very easily.




Do not use for 6.1011.030.



6.1011.130

Spare filter for RP 2 Guard/3.5, 10 pieces

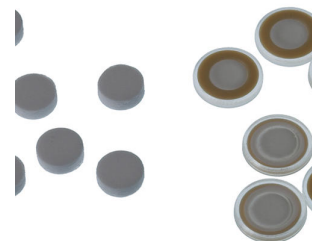

6.1031.500

Metrosep A Supp 16 Guard/4.0

The Metrosep A Supp 16 Guard/4.0 reliably protects the Metrosep A Supp 16 analytical separation columns against contamination. Thanks to the "On Column Guard System", the guard column is very easy to handle. The guard column screws easily onto the analytical column. No tools are required.



6.1031.510

Metrosep A Supp 16 S-Guard/4.0

The Metrosep A Supp 16 S-Guard/4.0 reliably protects the Metrosep A Supp 16 analytical separation columns against contamination. The guard column is connected to the separation column simply using capillary connections. No tools are required.

