



« ProfIC Vario 7 Anion » – système Professional IC Vario avec dilution inline et dialyse inline

ProfIC Vario 7 Anion

Le système Professional IC Vario doté de la dilution inline de Metrohm, de la dialyse inline, de la suppression séquentielle et de la détection de conductivité permet de déterminer de manière entièrement automatique les anions ou les cations dont la concentration est très élevée et dont la matrice est problématique.

Domaines d'application typiques :

- Eaux de processus, de lavage et usées
- Échantillons contenant de l'huile
- Extraction inline pour analyses de biodiesel

Représentation schématique

Contenu de la livraison ProfIC Vario 7 Anion

Qt.	Order no.	Description
-----	-----------	-------------


1 PCS

2.741.0010

741 Magnetic Stirrer

Magnetic stirrer for sample changers.

The magnetic stirrer is located beneath the sample rack.

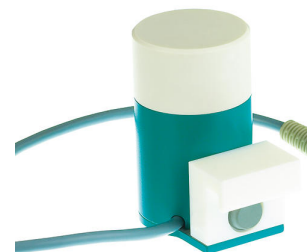

1 PCS

2.800.0010

800 Dosino

The 800 Dosino is a drive with write/read hardware for intelligent Dosing Units. With fixed cable (length 150 cm).



1 PCS

2.850.9010

IC Conductivity Detector

Compact and intelligent high performance conductivity detector for intelligent IC instruments. Outstanding temperature stability, the complete signal processing within the protected detector block and the latest generation of DSP – Digital Signal Processing – guarantee the highest precision of the measurement. No change of measuring ranges (not even automatic ones) is required, due to the dynamic working range.



1 PCS

2.858.0020

858 Professional Sample Processor – Pump

The 858 Professional Sample Processor – Pump processes samples from 500 µL to 500 mL. The sample transfer takes place either with the installed bidirectional two-channel peristaltic pump or with an 800 Dosino.



The 940 Professional IC Vario ONE/SeS/PP/Prep 1 is the intelligent IC instrument with sequential suppression and peristaltic pump for suppressor regeneration in combination with Metrohm Inline Sample Preparation, e.g., Inline Ultrafiltration or Inline Dialysis. The instrument can be used with any separation and detection methods.



Typical areas of application:

- Anion or cation determinations with sequential suppression after Inline Ultrafiltration or Inline Dialysis

Accessory set for Inline Dilution. For use with the 858 Professional Sample Processor, 800 Dosino and 741 Magnetic Stirrer.



Accessory set for rapid Inline Dialysis. For use with the 858 Professional Sample Processor and an additional two-channel peristaltic pump.



Accessoires optionnels

Order no.	Description	
2.941.0010	941 Eluent Production Module	
	The 941 Eluent Production Module enables the automatic production of an eluent. It enables continuous working without manual intervention and guarantees stable retention times. It can be combined with all Metrohm IC instruments and is monitored by MagIC Net.	
6.2041.440	Sample rack 148 x 11 mL + 3 x 300 mL	
	Sample rack for 148 x 11 mL samples with 3 rinsing beakers.	
6.2041.760	Sample rack 54 x 11 mL + 1 x 300 mL	
	Sample rack with small diameter for 54 x 11 mL sample vessels in two rows and additionally 1 x 300 mL bottle.	
6.2743.050	Sample vessel 11 mL	
	2000 pieces. For IC Sample Processors and VA Autosamplers. The IC sample vessels are made of polypropylene (PP) plastic for health-care applications. Production and packaging both take place under sterile and clean room conditions (Clean Room Class 7, in accordance with DIN ISO 14644). This guarantees minimal degrees of contamination. The sample vessels are subjected to random checks for leachable anions and cations.	

6.2743.070 Stopper with perforation
2000 pieces. For IC Sample Processors



6.2744.080 M6 thread / UNF 10/32 adapter
With M6 outer thread. For IC instruments.



6.2832.000 MSM Rotor A
Suppressor rotor for all IC instruments with MSM (Metrohm Suppressor Module)



6.2842.000 MSM-HC Rotor A
Suppressor rotor for all IC instruments with MSM-HC (Metrohm Suppressor Module with high capacity)



6.2842.020 Adapter sleeve, MSM

Not only the MSM Rotor A but also the MSM-LC Rotor A can be installed in the rotor housing with the aid of the adapter sleeve.



6.2842.200 MSM-HC Rotor C

Suppressor rotor for all IC instruments with MSM-HC drive (Metrohm Suppressor Module with high capacity) for cation suppression.



6.2844.000 MSM-LC Rotor A

Suppressor rotor for use with 2 mm columns.



6.3032.150 807 Dosing Unit 5 mL

807 Dosing Unit with integrated data chip with 5 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



Professional chromatography software from Metrohm for controlling all intelligent Professional IC systems, Compact IC systems and their peripherals, all detectors and the various Auto samplers, 800 Dosino, 771 Compact Interface, etc. The software permits checks, data acquisition, data evaluation and data monitoring as well as report generation of ion chromatographic analyses.



Graphical user interface for routine operations, extensive database programs, method development, configuration and manual control panel; very flexible user management, efficient database operations, extensive data export functions, individually configurable report generator, control and monitoring of all system components and the chromatography results.

MagIC Net Professional complies fully with FDA Regulation 21 CFR Part 11 as well as GLP.

MagIC Net is available in 15 dialog languages: German, English, Chinese, Traditional Chinese, French, Italian, Spanish, Portuguese, Bulgarian, Czech, Hungarian, Japanese, Korean, Russian, Slovakian

With this version, the graphic user interface has been modernized, which improves, among other things, the scaling and sharpness of the user interface on 4K monitors.

1 license

The installation and documents are supplied on a USB flash drive.

Client/server-compatible chromatography software from Metrohm for controlling all intelligent Professional IC systems, Compact IC systems and their peripherals, all detectors and the various Autosamplers, 800 Dosino, 771 Compact Interface, etc.

The software permits checks, data acquisition, data evaluation and data monitoring as well as report generation of ion chromatographic analyses. Data and method are centrally stored on a Server.



Graphical user interface for routine operations, extensive database programs, method development, configuration and manual control panel; very flexible user management, efficient database operations, extensive data export functions, individually configurable report generator, control and monitoring of all system components and the chromatography results.

MagIC Net Multi complies fully with FDA Regulation 21 CFR Part 11 as well as with GLP.

MagIC Net is available in 15 dialog languages: German, English, Chinese, Traditional Chinese, French, Italian, Spanish, Portuguese, Bulgarian, Czech, Hungarian, Japanese, Korean, Russian, Slovakian

With this version, the graphic user interface has been modernized, which improves, among other things, the scaling and sharpness of the user interface on 4K monitors.

Performance for client/server systems with servers outside the in-house company network has been improved.

Client/server version with 3 licenses
