


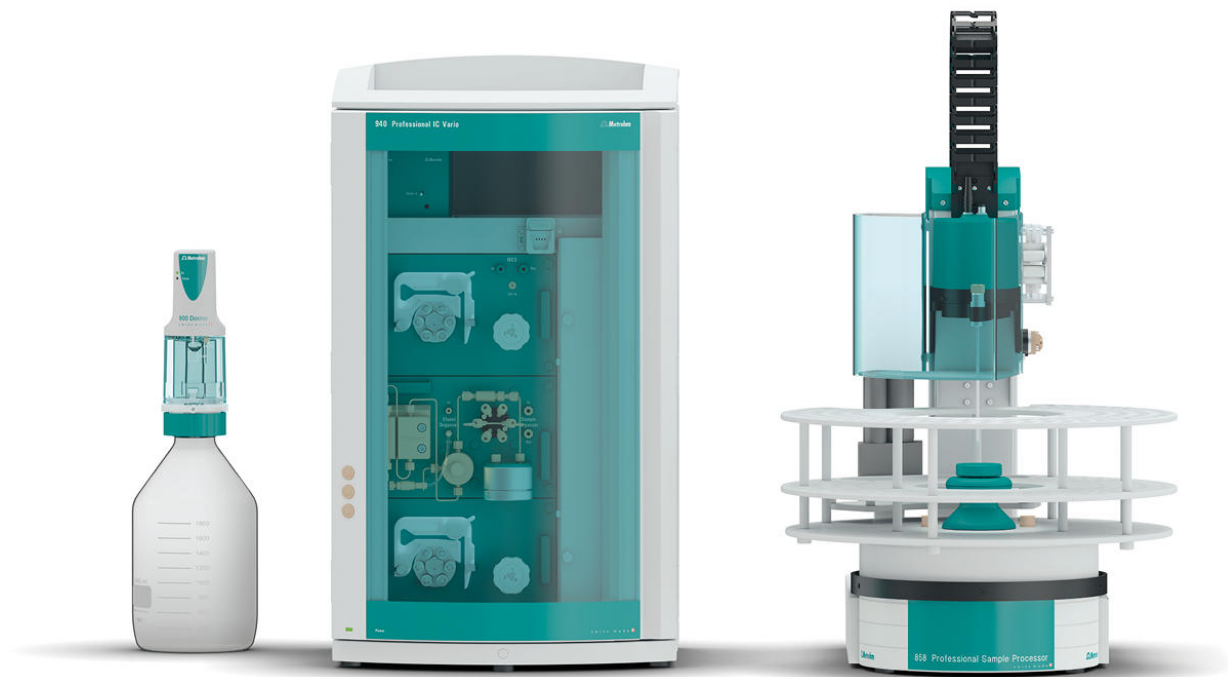

« ProfIC Vario 11 Anion » – système Professional IC Vario avec neutralisation inline, technique de préconcentration intelligente et élimination de la matrice inline

ProfIC Vario 11 Anion


Le système Professional IC Vario doté de la neutralisation inline de Metrohm, de la technique de préconcentration inline intelligente de Metrohm, de l'élimination de la matrice inline, de la suppression séquentielle et de la détection de conductivité permet d'effectuer des analyses des anions dans des matrices complexes à l'état d'ultra-traces jusque dans la gamme mg/L.

Domaines d'application typiques :

- Analyses de traces dans des bases et acides concentrés
- Analyses des ultra-traces pour des centrales nucléaires dans des matrices telles que hydroxyde de lithium/acide borique

Représentation schématique

Contenu de la livraison ProfIC Vario 11 Anion

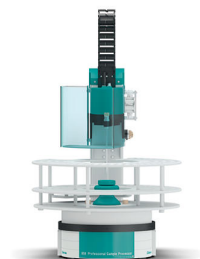Qt.	Order no.	Description
1 PCS	2.800.0010	800 Dosino
		The 800 Dosino is a drive with write/read hardware for intelligent Dosing Units. With fixed cable (length 150 cm).
		
1 PCS	2.850.9010	IC Conductivity Detector
		Compact and intelligent high performance conductivity detector for intelligent IC instruments. Outstanding temperature stability, the complete signal processing within the protected detector block and the latest generation of DSP – Digital Signal Processing – guarantee the highest precision of the measurement. No change of measuring ranges (not even automatic ones) is required, due to the dynamic working range.
		
1 PCS	2.858.0030	858 Professional Sample Processor – Pump – Injector
		The 858 Professional Sample Processor – Pump – Injector processes samples from 500 µL to 500 mL. The sample transfer takes place either with the installed bidirectional two-channel peristaltic pump or by means of an 800 Dosino. In addition, a six-way injection valve is available for the sample preparation.
		

The 940 Professional IC Vario ONE/SeS/PP/Prep 3 is the intelligent IC instrument with sequential suppression with a peristaltic pump for suppressor regeneration in combination with Metrohm Inline Sample Preparation, e.g., Inline Neutralization or Inline Cation Removal. The instrument can be used with any separation and detection methods.



Typical areas of application:

- Anion determinations with sequential suppression after Inline Neutralization
- Online trace analysis with Inline Neutralization and subsequent MiPCT-ME
- Anion analysis after cation removal in samples with high transition metal concentrations, e.g., electroplating baths

The Metrosep I Trap 1 - 100/4.0 column is used to eliminate ionic, i.e. cationic and anionic, contaminants from aqueous solutions. Its use is especially recommended for the purification of the transfer water in combination with "MISP" or "Metrohm Inline Sample Preparation". Using the Metrosep I Trap 1 - 100/4.0 column can significantly reduce the influence of the transfer water on the system blank.



Outer thread M8, inner thread M6, e.g. for mounting 6.2418.XXX Air tubes in 6.1429.XXX Reaction vessels of stability measuring instruments.



1 PCS

6.2832.000

MSM Rotor A

Suppressor rotor for all IC instruments with MSM (Metrohm Suppressor Module)



1 PCS

6.2842.020

Adapter sleeve, MSM

Not only the MSM Rotor A but also the MSM-LC Rotor A can be installed in the rotor housing with the aid of the adapter sleeve.



1 PCS

6.5330.140

IC equipment: MiPCT

Accessory set for Inline Preconcentration (MiPCT) including Liquid Handling Station, left.



Accessoires optionnels

Order no.	Description	
2.941.0010	941 Eluent Production Module	
	The 941 Eluent Production Module enables the automatic production of an eluent. It enables continuous working without manual intervention and guarantees stable retention times. It can be combined with all Metrohm IC instruments and is monitored by MagIC Net.	
6.1608.100	PE bottle 50 mL	
	Bottle for Extension Module Liquid Handling.	
6.2744.080	M6 thread / UNF 10/32 adapter	
	With M6 outer thread. For IC instruments.	
6.2820.000	Al foil	
	1000 circular sheets of 80 mm diameter. Thickness: 0.010 mm. Used together with 6.2037.030 Foil holder for 6.1432.210 Sample beaker. For 674 KF Sample Changer	

6.2842.000

MSM-HC Rotor A

Suppressor rotor for all IC instruments with MSM-HC (Metrohm Suppressor Module with high capacity)



6.2842.200

MSM-HC Rotor C

Suppressor rotor for all IC instruments with MSM-HC drive (Metrohm Suppressor Module with high capacity) for cation suppression.



6.2844.000

MSM-LC Rotor A

Suppressor rotor for use with 2 mm columns.



6.3032.150

807 Dosing Unit 5 mL

807 Dosing Unit with integrated data chip with 5 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



Professional chromatography software from Metrohm for controlling all intelligent Professional IC systems, Compact IC systems and their peripherals, all detectors and the various Auto samplers, 800 Dosino, 771 Compact Interface, etc. The software permits checks, data acquisition, data evaluation and data monitoring as well as report generation of ion chromatographic analyses.



Graphical user interface for routine operations, extensive database programs, method development, configuration and manual control panel; very flexible user management, efficient database operations, extensive data export functions, individually configurable report generator, control and monitoring of all system components and the chromatography results.

MagIC Net Professional complies fully with FDA Regulation 21 CFR Part 11 as well as GLP.

MagIC Net is available in 15 dialog languages: German, English, Chinese, Traditional Chinese, French, Italian, Spanish, Portuguese, Bulgarian, Czech, Hungarian, Japanese, Korean, Russian, Slovakian

With this version, the graphic user interface has been modernized, which improves, among other things, the scaling and sharpness of the user interface on 4K monitors.

1 license

The installation and documents are supplied on a USB flash drive.

Client/server-compatible chromatography software from Metrohm for controlling all intelligent Professional IC systems, Compact IC systems and their peripherals, all detectors and the various Autosamplers, 800 Dosino, 771 Compact Interface, etc.

The software permits checks, data acquisition, data evaluation and data monitoring as well as report generation of ion chromatographic analyses. Data and method are centrally stored on a Server.



Graphical user interface for routine operations, extensive database programs, method development, configuration and manual control panel; very flexible user management, efficient database operations, extensive data export functions, individually configurable report generator, control and monitoring of all system components and the chromatography results.

MagIC Net Multi complies fully with FDA Regulation 21 CFR Part 11 as well as with GLP.

MagIC Net is available in 15 dialog languages: German, English, Chinese, Traditional Chinese, French, Italian, Spanish, Portuguese, Bulgarian, Czech, Hungarian, Japanese, Korean, Russian, Slovakian

With this version, the graphic user interface has been modernized, which improves, among other things, the scaling and sharpness of the user interface on 4K monitors.

Performance for client/server systems with servers outside the in-house company network has been improved.

Client/server version with 3 licenses

Sample rack for 858 Professional Sample Processor, 815 Robotic USB Sample Processor XL, suitable for 35 sample vessels and 3 rinsing bottles. The following sample vessels and rinsing bottles can be used: 6.1608.100, 60 mL and 250 mL Nalgene FEP vessels.

Plastic: Polyvinyl chloride (PVC)

