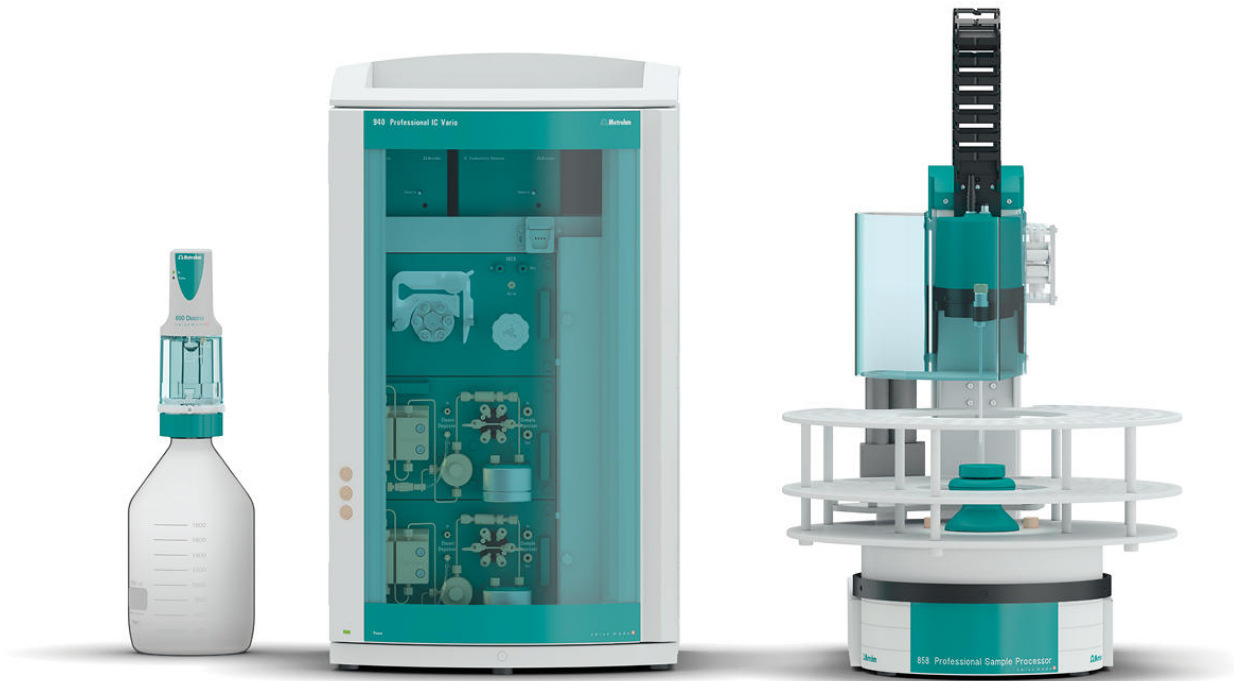




# « ProfIC Vario 14 AnCat » – système Professional IC Vario avec technique intelligente d'injection pick-up de Metrohm (MiPuT)

ProfIC Vario 14 AnCat

Le système Professional IC Vario doté de la technique intelligente d'injection pick-up de Metrohm (MiPuT), de la suppression séquentielle et de la détection de conductivité permet de déterminer des anions et des cations à partir de volumes d'échantillon extrêmement faibles.

Domaines d'application typiques :

- Utilisation universelle pour tous les types d'échantillons, si seuls de faibles volumes d'échantillon sont disponibles
- Analyse de processus biochimiques
- Extraits d'échantillons d'environnement

Représentation schématique

Contenu de la livraison ProfIC Vario 14 AnCat

| Qt. | Order no. | Description |
|-----|-----------|-------------|
|-----|-----------|-------------|

1 PCS

2.800.0010

800 Dosino

The 800 Dosino is a drive with write/read hardware for intelligent Dosing Units. With fixed cable (length 150 cm).



2 PCS

2.850.9010

IC Conductivity Detector

Compact and intelligent high performance conductivity detector for intelligent IC instruments. Outstanding temperature stability, the complete signal processing within the protected detector block and the latest generation of DSP – Digital Signal Processing – guarantee the highest precision of the measurement. No change of measuring ranges (not even automatic ones) is required, due to the dynamic working range.



1 PCS

2.858.0020

858 Professional Sample Processor – Pump

The 858 Professional Sample Processor – Pump processes samples from 500 µL to 500 mL. The sample transfer takes place either with the installed bidirectional two-channel peristaltic pump or with an 800 Dosino.



1 PCS

2.940.2500

940 Professional IC Vario TWO/SeS/PP

The 940 Professional IC Vario TWO/SeS/PP is the intelligent two-channel IC instrument with sequential suppression (one channel) and a peristaltic pump for suppressor regeneration. The instrument can be used with any separation and detection methods.

Typical areas of application:

- Standard instrument for parallel anion and cation determinations
- Trace analysis for anions and cations
- Online monitoring for anions and cations



1 PCS

6.1825.290

PEEK sample loop 250 µL

Sample loop for «Metrohm intelligent Partial loop Technique» (MiPT) and for PILS (particle into liquid sampler). With 2 Peek pressure screws.



---

1 PCS

6.5330.170

IC equipment: MiPuT

Accessory set for assembling a Dosino for sample pick-up (Metrohm intelligent Pick-up Injection Technique).



## Accessoires optionnels

| Order no.  | Description                                  |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                   |
|------------|----------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| 2.941.0010 | 941 Eluent Production Module                 | <p>The 941 Eluent Production Module enables the automatic production of an eluent. It enables continuous working without manual intervention and guarantees stable retention times. It can be combined with all Metrohm IC instruments and is monitored by MagIC Net.</p>                                                                                                                                                                                                                                      |
| 6.1050.420 | Metrosep C 4 - 150/4.0                       | <p>The Metrosep C 4 - 150/4.0 is the universal standard column in cation analysis for accomplishing high separating efficiency in a short time. The Metrosep C 4 - 150/4.0 is the ideal separation column for the analysis of alkaline and earth alkaline metals in aqueous media.</p>                                                                                                                                                                                                                         |
| 6.1050.500 | Metrosep C 4 Guard/4.0                       | <p>The Metrosep C 4 Guard/4.0 contains the C-4 column material and is used to protect all Metrosep cation columns that have a substrate based on silica gel. Particles and contamination are reliably retained, considerably prolonging the service life of the analytical separation column. The economical price is an additional plus. The Metrosep C 4 Guard/4.0 also works based on the "On Column Guard System" and is screwed directly onto the separation column with virtually no dead volume.</p>  |
| 6.2041.480 | Sample rack, 159 x 2 mL and 3 x 300 mL (PVC) | <p>Sample rack for 858 Professional Sample Processor, suitable for 159 sample vessels and 3 rinsing bottles. The following sample vessels can be used: 2 mL (Eppendorf), 6.1608.080, 6.1627.000.</p> <p>Plastic: Polyvinyl chloride (PVC)</p>                                                                                                                                                                                                                                                                |

6.2769.110 EMPTY sensor for level sensor for 2 L bottles

This sensor, mounted in the 6.1626.000 bottle attachment, results in the level sensor EMPTY for 2 L bottles. Included are two sensor rods, one aspiration tubing and two holders for fixing the tubing to one of the sensor rods.



---

6.2832.000 MSM Rotor A

Suppressor rotor for all IC instruments with MSM (Metrohm Suppressor Module)



---

6.2842.000 MSM-HC Rotor A

Suppressor rotor for all IC instruments with MSM-HC (Metrohm Suppressor Module with high capacity)



---

6.2842.020 Adapter sleeve, MSM

Not only the MSM Rotor A but also the MSM-LC Rotor A can be installed in the rotor housing with the aid of the adapter sleeve.



6.2844.000

MSM-LC Rotor A

Suppressor rotor for use with 2 mm columns.



---

6.5330.090

IC equipment: Additional eluent on the Eluent Production Module

Accessory set for Eluent Production Module. Contains all parts necessary for expansion by one eluent. The level sensor must be ordered in accordance with the application.



Professional chromatography software from Metrohm for controlling all intelligent Professional IC systems, Compact IC systems and their peripherals, all detectors and the various Auto samplers, 800 Dosino, 771 Compact Interface, etc. The software permits checks, data acquisition, data evaluation and data monitoring as well as report generation of ion chromatographic analyses.



Graphical user interface for routine operations, extensive database programs, method development, configuration and manual control panel; very flexible user management, efficient database operations, extensive data export functions, individually configurable report generator, control and monitoring of all system components and the chromatography results.

MagIC Net Professional complies fully with FDA Regulation 21 CFR Part 11 as well as GLP.

MagIC Net is available in 15 dialog languages: German, English, Chinese, Traditional Chinese, French, Italian, Spanish, Portuguese, Bulgarian, Czech, Hungarian, Japanese, Korean, Russian, Slovakian

With this version, the graphic user interface has been modernized, which improves, among other things, the scaling and sharpness of the user interface on 4K monitors.

1 license

The installation and documents are supplied on a USB flash drive.

Client/server-compatible chromatography software from Metrohm for controlling all intelligent Professional IC systems, Compact IC systems and their peripherals, all detectors and the various Autosamplers, 800 Dosino, 771 Compact Interface, etc.

The software permits checks, data acquisition, data evaluation and data monitoring as well as report generation of ion chromatographic analyses. Data and method are centrally stored on a Server.



Graphical user interface for routine operations, extensive database programs, method development, configuration and manual control panel; very flexible user management, efficient database operations, extensive data export functions, individually configurable report generator, control and monitoring of all system components and the chromatography results.

MagIC Net Multi complies fully with FDA Regulation 21 CFR Part 11 as well as with GLP.

MagIC Net is available in 15 dialog languages: German, English, Chinese, Traditional Chinese, French, Italian, Spanish, Portuguese, Bulgarian, Czech, Hungarian, Japanese, Korean, Russian, Slovakian

With this version, the graphic user interface has been modernized, which improves, among other things, the scaling and sharpness of the user interface on 4K monitors.

Performance for client/server systems with servers outside the in-house company network has been improved.

Client/server version with 3 licenses

---