



907 Titrando

2.907.0010

Titolatore di fascia alta per la titolazione Karl Fischer potenziometrica e volumetrica con un'interfaccia di misura e unità di dosaggio Dosino.

- fino a quattro sistemi di dosaggio del tipo 800 Dosino
- titolazione dinamica (DET), monotonica (MET) e a punto finale (SET), titolazioni enzimatiche e pH-STAT (STAT), titolazione Karl Fischer (KFT)
- Misura con elettrodi iono-selettivi (MEAS CONC)
- elettrodi intelligenti "iTrode"
- Funzioni di dosaggio con monitoraggio, trasferimento liquidi
- quattro connettori MSB per ulteriori agitatori o sistemi di dosaggio
- connettore USB
- Utilizzo con software OMNIS, software *tiamo* o Touch Control
- Conforme ai requisiti GMP/GLP e FDA, nonché 21 CFR Parte 11, se necessario

Parti incluse 2.907.0010

Qt.	Order no.	Descrizione
-----	-----------	-------------

High-end titrator for use with intelligent electrodes – iTrodes – with


- up to four dosing systems of the 800 Dosino type
- dynamic (DET), monotonic (MET) and endpoint titration (SET), enzymatic and pH-STAT titrations (STAT), Karl Fischer titration (KFT)
- Measurement with ion-selective electrodes (MEAS CONC)
- Dosing functions with monitoring (DOS), Liquid Handling
- four MSB connectors
- one galvanically isolated measuring interface
- USB connector



Holding spring for reagent bottles in exchange units.



Accessori opzionali

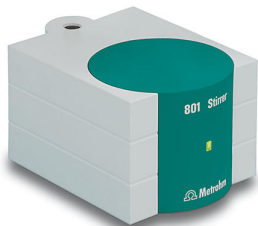Order no.	Descrizione	
2.800.0010	800 Dosino	The 800 Dosino is a drive with write/read hardware for intelligent Dosing Units. With fixed cable (length 150 cm). 
2.801.0040	801 Stirrer with stand	Magnetic stirrer including base plate, support rod and electrode holder for supplementing the Titrino plus, Dosimat plus, Titrandos, Sample Processors, 805 Dosimat and 780/781 pH meters as well as the 856 and 867 measuring modules. With permanently attached cable for MSB (Metrohm Serial Bus). 
2.802.0040	802 Stirrer for 804 Ti Stand	Rod stirrer including 6.1909.010 stirring propeller. 
2.804.0040	804 Ti Stand with stand	Titration stand and controller for 802 Rod Stirrer. The 804 Ti Stand together with the optional 802 Rod Stirrer provides an alternative to the magnetic stirrer. Ti Stand with base plate, support rod and electrode holder. 

2.805.0010

805 Dosimat

Dosing device for the Titrando and Sample Processors with read/write hardware for intelligent Exchange Units. With permanently attached cable. Without Exchange Unit.



2.846.0010

846 Dosing Interface

USB-capable control unit for connecting a maximum of four 800 Dosinos or 805 Dosimats for dosing and Liquid Handling tasks. A Touch Control or the connection to a PC with OMNIS Software, tiamo™, MagIC Net, viva, or 797 VA Computrace is required for operation.



2.854.0010

854 iConnect


The 854 iConnect – Electrode cable and measuring amplifier for intelligent "iTrodes" electrodes (cable length 150 cm).



2.900.0010

900 Touch Control

Operating unit for the Titrandos, USB Sample Processors, 856 Conductivity Module, 867 pH Module and 846 Dosing Interface. Touch-sensitive, high-resolution color display, simple and intuitive operation, thanks to Favorites for direct method access. With integrated Ethernet interface for direct connection to the Internet and USB interface for connecting USB printers or a USB memory stick.

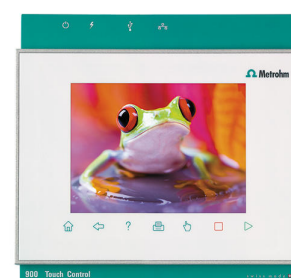

Dialog languages: German, English, Chinese, French, Spanish, Portuguese, Russian, Korean, Polish and Italian.

6.2103.130 Adapter red 2 mm plug / 4 mm socket
For connecting plug B (4 mm) to 2 mm socket.



6.2103.140 Adapter black 2 mm plug / B socket 4 mm
For connecting plug B (4 mm) to 2 mm socket.



6.2104.020 Electrode cable / 1 m / F
For connecting electrodes with Metrohm plug-in head G to Metrohm instruments (socket F).



6.2104.600 Electrode cable for plug in head U/plug F, 2x2 mm B, 1m
For connecting electrodes with Metrohm plug-in head U to Metrohm instruments (socket F).



6.2148.010

Remote Box MSB

Additional remote interface for the connection of devices that can be controlled via remote lines. With permanently attached cable.



6.2148.020

RS-232/USB box

Extends the Titrando, USB Sample Processors, Dosing Interface, Titrosampler by two RS-232 interfaces for connecting e.g. balances, PC, LIMS. Requires 6.2151.030 cable.



6.2168.000

Mini USB adapter cable, 25 cm

Mini USB adapter cable for connecting one 854 iConnect to an 867 pH Module, 901 to 907 Titrandos and 916 Ti-Touch, 25 cm. The adapter cable remains plugged into the instrument as strain relief and the 854 iConnect is only plugged into and unplugged from the adapter cable.



6.3026.110

Exchange Unit 1 mL

Exchange Unit with integrated data chip with 1 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, antidiffusion buret tip and standard amber glass reagent bottle.



6.3026.150

Exchange Unit 5 mL

Exchange Unit with integrated data chip with 5 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, anti-diffusion buret tip and standard amber glass reagent bottle.



6.3026.210

Exchange Unit 10 mL

Exchange unit with integrated data chip with 10 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, anti-diffusion buret tip and standard amber glass reagent bottle.



6.3026.220

Exchange Unit 20 mL

Exchange unit with integrated data chip with 20 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, anti-diffusion buret tip and standard amber glass reagent bottle.



6.3026.250

Exchange unit 50 mL

Exchange unit with integrated data chip with 50 mL glass cylinder and light protection. PCTFE/PTFE flat stopcock, FEP tubing connection, antidiffusion tip and standard amber glass reagent bottle.



6.3032.120

807 Dosing Unit 2 mL

807 Dosing Unit with integrated data chip with 2 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.3032.150

807 Dosing Unit 5 mL

807 Dosing Unit with integrated data chip with 5 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.3032.210

807 Dosing Unit 10 mL

807 Dosing Unit with integrated data chip with 10 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.3032.220

807 Dosing Unit 20 mL

807 Dosing Unit with integrated data chip with 20 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.



6.3032.250

807 Dosing Unit 50 mL

807 Dosing Unit with integrated data chip with 50 mL glass cylinder and light protection, mountable to a reagent bottle with ISO/DIN GL 45 glass thread. FEP tubing connection, antidiffusion tip.

