



913 pH/DO Meter

2.913.0020

913 pH/DO Meter pH2

- 2
- pH pH
- O2 LumitrodepH O2
- (IP67)
- LCD
- USB
- (10,000)
- PIN
- IDGLP

2.913.0020

Qt.	Order no.
-----	-----------

The 913 pH/DO Meter offers two measuring inputs for determining the pH value and dissolved oxygen in water.

- Portable meter with built-in battery pack and two galvanically isolated measuring inputs.
- Analog pH measuring input for standard pH electrodes
- Digital measuring input for O₂ measurement using O₂ Lumitrode or intelligent pH electrodes
- Robust, water-tight, and dust-tight housing (IP67) for tough outdoor and laboratory use
- LCD color display with background illumination making results easy to read
- USB interface for simple data export to PC or printer
- Large internal memory (10,000 data sets)
- Pin-protected User mode and Expert mode, prevents unwanted parameter changes
- GLP-compliant printout and data export with User ID and timestamp



Practical holder for fastening the electrode storage vessels to the 912, 913 or 914 pH/Conductometers.



Carrying strap for 912/913/914 meters



1 PCS

6.2151.100

Adapter USB MINI (OTG) - USB A

For connecting USB instruments.



1 PCS

6.2151.110

Metrohm USB Mini B cable (OTG) - USB A, 1.8 m

For connecting USB instruments.



1 PCS

6.2166.100

USB power supply unit 5.25 V / 1.53 A

USB power supply unit for 912 / 913 / 914

Efficiency Level VI



Order no.**2.142.0100 Custom Q3X thermal printer**

Compact printer with USB interface for

- 900 Touch Control
- 915 KF Ti-Touch
- 916 Ti-Touch
- 917 Coulometer
- 877 / 848 Titrino plus
- 865 / 876 Dosimat plus
- 91X Meter (cable 6.2151.140)
- Eco Dosimat / Titrator
- 862 Compact Titrosampler
- 870 KF Titrino plus
- 899 Coulometer



Paper width 60 mm (40 characters). Including 6.2151.120 USB cable.

2.854.0010 854 iConnect

The 854 iConnect – Electrode cable and measuring amplifier for intelligent "iTrodes" electrodes (cable length 150 cm).



6.00226.600 Spearhead electrode with Pt1000

Maintenance-free combined pH electrode (gel electrolyte) for piercing measurements of all types (e.g., with cheese, meat, dough) with integrated Pt1000 temperature sensor. The electrode is stored in saturated potassium chloride solution $c(\text{KCl}) = \text{sat.}$ (6.2308.000) and is not suitable for low-ion solutions. The ageing indicator gives early indication of when the electrode needs to be replaced.



Combined pH electrode for measurements in very small sample volumes (>50 µL) and biological samples.

Idrolyte (6.2308.040) is used as reference electrolyte and storage solution.



Combined pH electrode with integrated Pt1000 temperature sensor and fixed cable (1.2 m). This electrode is suitable for routine pH measurements in solutions that do not contain precipitates, proteins, or sulfides. This electrode is mechanically resistant thanks to the robust /unbreakable plastic shaft made of polypropylene and impact protection for the glass membrane.

Reference electrolyte: c(KCl) = 3 mol/L, storage in storage solution.



Combined pH electrode with integrated temperature sensor (Pt1000) and fixed cable conforming to IP67 (1.2 m).

This electrode is suitable for an introduction to GLP-compliant pH measurements in solutions that do not contain precipitates, proteins, or sulfides. This electrode is mechanically resistant thanks to the robust /unbreakable plastic shaft made of polypropylene and impact protection for the glass membrane.

Additionally, this electrode offers a waterproof plug in accordance with IP67 for mobile use with Metrohm pH meters.

Reference electrolyte: c(KCl) = 3 mol/L, storage in storage solution.




Combined pH electrode with integrated Pt1000 temperature sensor and fixed cable (1.2 m, diameter banana plug 2 mm). This electrode is particularly suitable

- for pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- for long-term measurements

The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

Alternatively: reference electrolyte for measurements at $T > 80^\circ\text{C}$: Idrolyte, storage in Idrolyte.

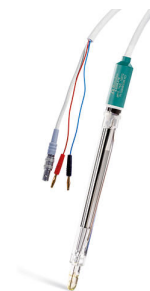

Intelligent, combined pH electrode with integrated memory chip for storing sensor data and Pt1000 temperature sensor. This electrode is particularly suitable:

- for pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- for long-term measurements

The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

Alternatively: reference electrolyte for measurements at $T > 80^\circ\text{C}$: Idrolyte, storage in Idrolyte.

iTrodes can be connected to Titrando, Ti-Touch or 913/914 meters.



The optical sensor for measuring dissolved oxygen (DO) can be used with a 913 pH/DO meter or with a 914 pH/DO conductometer. The measuring principle of the sensor is based on luminescence quenching. The space-saving and maintenance-free sensor is suitable for DO measurement, for example in:

- Water quality control
- Wastewater industry
- Beverage production
- Fish farming

This sensor is supplied with a calibration vessel.

Where necessary, it is easy to replace the measurement cap (O₂ cap) that contains the oxygen-intensive luminophore.



Stand plate for converting a mobile 912/913/914 pH/Conductometer into a laboratory meter.



For connecting electrodes with Metrohm plug-in head U to Metrohm instruments (socket F).



6.2151.140

Y cable USB A St - USB B St - Mini B St

Y cable for connection of a USB printer to the pH/Conductometers 912 / 913 / 914. This cable allows a printer and the power supply unit to be connected to the measuring instrument at the same time.



6.2166.500

12 V USB adapter for 912 / 913 / 914 pH/Conductometer

12 V USB adapter for 912 / 913 / 914 pH/Conductometer.



6.2307.230

Buffer solutions pH 4, 7 and 9

Mixed buffer solutions pH 4.00/7.00/9.00 (25 °C) in single use sachets, colourless, box of 3 x 10 x 30 mL



6.2313.000

Electrolyte 3 mol/L KCl (1000 mL)

Electrolyte solution $c(\text{KCl}) = 3 \text{ mol/L}$ (for Ag/AgCl reference systems)



6.2325.000

pHit kit

Maintenance kit for pH electrodes

The kit contains:

- 50 mL cleaning solution
- 50 mL 3M KCl solution
- 50 mL storage solution
- 2 Storage vessels
- Instructions for use



6.2329.000

Oxygen standard 0%

Oxygen standard for 0% calibration of the O₂-Lumitrode. Contains 5 sachets of 30 mL each.



6.5623.000

O₂ Lumitrode replacement set

Replacement set for the O₂ Lumitrode. The set contains a certified O₂ cap as replacement, a calibration vessel, and 5 x 30 mL oxygen standard, 0%.

