



913 pH Meter

2.913.0110

pH/mV 2pH

- pH2pH
- MetrohmpH pH
- MetrohmpH pH
- (IP67)
- LCD
- USB
- (10000)
- PIN
- IDGLP
-

2.913.0110

Qt.

Order no.

pH/mV and temperature measurement for routine use in the laboratory and on the road.

pH or mV and temperature can be measured in parallel with one analog and one intelligent sensor and output to a large colour display. Important information such as charge state, user, IDs can be clearly seen at a glance. A PIN-protected expert mode protects against unintentional changes of different parameters.

The meter is furnished with an accumulator for mobile use that can be charged practically anywhere. It naturally also satisfies the requirements of IP67.

A stand plate allows the mobile meter to be easily converted into a laboratory meter and vice versa.

Very large measured value memory (10,000 data sets) and USB interface (GLP-compliant printout or data export with optional management of the data in tiBase) offer professional data handling.



Combined pH electrode with integrated temperature sensor (Pt1000) and fixed cable conforming to IP67 (1.2 m).

This electrode is suitable for an introduction to GLP-compliant pH measurements in solutions that do not contain precipitates, proteins, or sulfides. This electrode is mechanically resistant thanks to the robust /unbreakable plastic shaft made of polypropylene and impact protection for the glass membrane.

Additionally, this electrode offers a waterproof plug in accordance with IP67 for mobile use with Metrohm pH meters.

Reference electrolyte: c(KCl) = 3 mol/L, storage in storage solution.



2 PCS

6.1613.010

Bottle / 25 mL



1 PCS

6.1613.020

Bottle / 25 mL / label pH 7 (empty)



1 PCS

6.1613.030

Bottle / 25 mL / label pH 4 (empty)



1 PCS

6.1614.000

Wash bottle / 250 mL



1 PCS

6.2008.060

Holder for electrode storage vessel

Practical holder for fastening the electrode storage vessels to the 912, 913 or 914 pH/Conductometers.



1 PCS

6.2050.010

Carrying strap for 912/913/914

Carrying strap for 912/913/914 meters



1 PCS

6.2151.100

Adapter USB MINI (OTG) - USB A

For connecting USB instruments.



1 PCS

6.2151.110

Metrohm USB Mini B cable (OTG) - USB A, 1.8 m

For connecting USB instruments.



1 PCS

6.2166.100

USB power supply unit 5.25 V / 1.53 A

USB power supply unit for 912 / 913 / 914

Efficiency Level VI



1 PCS

6.2307.230

Buffer solutions pH 4, 7 and 9

Mixed buffer solutions pH 4.00/7.00/9.00 (25 °C) in single use sachets, colourless, box of 3 x 10 x 30 mL



1 PCS

6.2308.050

Electrolyte 3 mol/L KCl (50 mL)

Electrolyte solution $c(\text{KCl}) = 3 \text{ mol/L}$ (for Ag/AgCl reference systems)



1 PCS

6.2716.060

Case for 912 / 913 / 914

Case for the 912 / 913 / 914 pH/Conductometers.



1 PCS

6.2717.000

PP beaker, 100 mL

Beaker made of PP, 100 mL.



Order no.**2.142.0100 Custom Q3X thermal printer**

Compact printer with USB interface for

- 900 Touch Control
- 915 KF Ti-Touch
- 916 Ti-Touch
- 917 Coulometer
- 877 / 848 Titrino plus
- 865 / 876 Dosimat plus
- 91X Meter (cable 6.2151.140)
- Eco Dosimat / Titrator
- 862 Compact Titrosampler
- 870 KF Titrino plus
- 899 Coulometer



Paper width 60 mm (40 characters). Including 6.2151.120 USB cable.

2.854.0010 854 iConnect

The 854 iConnect – Electrode cable and measuring amplifier for intelligent "iTrodes" electrodes (cable length 150 cm).



6.00226.600 Spearhead electrode with Pt1000

Maintenance-free combined pH electrode (gel electrolyte) for piercing measurements of all types (e.g., with cheese, meat, dough) with integrated Pt1000 temperature sensor. The electrode is stored in saturated potassium chloride solution $c(\text{KCl}) = \text{sat.}$ (6.2308.000) and is not suitable for low-ion solutions. The ageing indicator gives early indication of when the electrode needs to be replaced.




Combined pH electrode for measurements in very small sample volumes (>50 µL) and biological samples.

Idrolyte (6.2308.040) is used as reference electrolyte and storage solution.



Combined pH electrode with integrated Pt1000 temperature sensor and fixed cable (1.2 m). This electrode is suitable for routine pH measurements in solutions that do not contain precipitates, proteins, or sulfides. This electrode is mechanically resistant thanks to the robust /unbreakable plastic shaft made of polypropylene and impact protection for the glass membrane.

Reference electrolyte: c(KCl) = 3 mol/L, storage in storage solution.

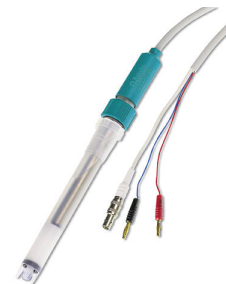


Combined pH electrode with integrated Pt1000 temperature sensor and fixed cable (1.2 m, diameter banana plug 2 mm). This electrode is particularly suitable

- for pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- for long-term measurements

The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: c(KCl) = 3 mol/L, storage in storage solution.

Alternatively: reference electrolyte for measurements at $T > 80^{\circ}\text{C}$: Idrolyte, storage in Idrolyte.

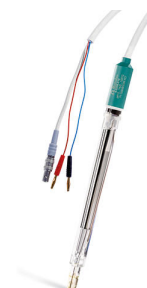

Intelligent, combined pH electrode with integrated memory chip for storing sensor data and Pt1000 temperature sensor. This electrode is particularly suitable:

- for pH measurements and titrations in difficult, viscous, or alkaline samples
- at elevated temperatures
- for long-term measurements



The fixed ground-joint diaphragm is insensitive to contamination.

Reference electrolyte: $c(\text{KCl}) = 3 \text{ mol/L}$, storage in storage solution.

Alternatively: reference electrolyte for measurements at $T > 80^\circ\text{C}$: Idrolyte, storage in Idrolyte.

iTrodes can be connected to Titrando, Ti-Touch or 913/914 meters.

6.2001.130**Stand plate for 912/913/914**

Stand plate for converting a mobile 912/913/914 pH/Conductometer into a laboratory meter.



6.2104.600**Electrode cable for plug in head U/plug F, 2x2 mm B, 1m**

For connecting electrodes with Metrohm plug-in head U to Metrohm instruments (socket F).



6.2151.140

Y cable USB A St - USB B St - Mini B St

Y cable for connection of a USB printer to the pH/Conductometers 912 / 913 / 914. This cable allows a printer and the power supply unit to be connected to the measuring instrument at the same time.



6.2166.500

12 V USB adapter for 912 / 913 / 914 pH/Conductometer

12 V USB adapter for 912 / 913 / 914 pH/Conductometer.



6.2313.000

Electrolyte 3 mol/L KCl (1000 mL)

Electrolyte solution $c(\text{KCl}) = 3 \text{ mol/L}$ (for Ag/AgCl reference systems)



6.2325.000

pHit kit

Maintenance kit for pH electrodes

The kit contains:

- 50 mL cleaning solution
- 50 mL 3M KCl solution
- 50 mL storage solution
- 2 Storage vessels
- Instructions for use

