



## Ag Titrode

6.00430.100

Combined silver ring electrode with a pH glass membrane as reference electrode.

This maintenance-free electrode is suitable for precipitation titrations when the pH value remains constant (titrant silver nitrate), for example of:

- chloride, bromide, iodide
- sulfides
- hydrogen sulfide
- mercaptans
- cyanides

This electrode is stored in distilled water.

Depending on the application, we recommend using an Ag Titrode with a coating (AgBr, AgCl or Ag<sub>2</sub>S coating), which can be ordered accordingly.

Phạm vi giao hàng 6.00430.100

Qt.	Order no.	Description
-----	-----------	-------------

1 PCS

6.1236.050

Sleeve with SGJ 14/12 mm

Sleeve with SGJ 14/12 mm, polyethylene.



---

1 PCS

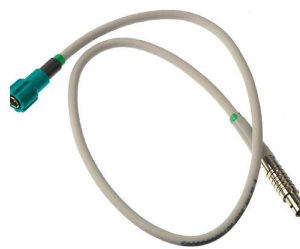

6.2008.040

Storage vessel

Together with 6.2008.050 storage vessel holder. it provides a support for the electrode on 807 Dosing Units.



## Phụ kiện tùy chọn

Order no.	Description	
6.02104.000	Electrode cable plug-in head G / plug P, 0.55 m  Cable for connecting an analog electrode (plug-in head G) to an analog measuring module from OMNIS.	
6.02104.010	Electrode cable plug-in head G / plug P, 1.5 m  Cable for connecting an analog electrode (plug-in head G) to an analog measuring module from OMNIS.	
6.2104.020	Electrode cable / 1 m / F  For connecting electrodes with Metrohm plug-in head G to Metrohm instruments (socket F).	